



KINEXUS HHD PROVIDER PORTAL USER MANUAL

© 2025 Fresenius USA, Inc. All Rights Reserved.

This document contains proprietary information of Fresenius Medical Care Renal Therapies Group, LLC and its affiliates (“Fresenius Medical Care”). The contents of this document may not be disclosed to third parties, copied, or duplicated in any form, in whole or in part, without the prior written permission of Fresenius Medical Care.

Fresenius Medical Care, the triangle logo, NxStage, Nx2me, Nx2me Connected Health, and kinexus are trademarks of Fresenius Medical Care Holdings, Inc. or its affiliated companies. All other trademarks are the property of their respective owners.

For more information, contact:

**Fresenius Medical Care AG & Co. KGaA
Else-Kröner-Str. 1
61352 Bad Homburg
Germany**

**Telephone: +49 6172 609-0
www.fmc-ag.com**

Table of Contents

Chapter 1: Indications for Use	6
Overview	6
Safety Notice	6
System Requirements	7
Screen Size	7
Operating System and Browser	7
Setup	8
Initial Onboarding and Customer Service	8
Customisation of the Portal	8
Support Resources	9
Customer Service	9
Technical Support	9
Warnings	9
Symbols and conventions	10
Accessing the User Manual	11
Logging into kinexus HHD Portal	13
Multi-Factor Authentication (MFA)	14
Password Reset	16
Steps to Reset Your Password	17
Chapter 2: Introduction to kinexus HHD	18
Overview	18
kinexus HHD Patient App (Nx2me)	19
kinexus HHD Provider Portal	19
Use	19
Disclaimer	19
User Roles	20
Chapter 3: Clinic Dashboard	21
Overview	21
Summary	22
Access <i>Clinic Overview Card</i>	23
All	23
Therapy Gaps	24

First 90 Days	24
Alerts & Alarms – Last 7 Days.....	25
Last 7 Days Treatments	26
Search by patient	26
Last 7 Days Treatments Table	26
How to Review Treatments	28
Treatments	29
Sort Treatment list.....	30
Patients	31
Sort Patients list	32
Patients table	32
+ Add Patient.....	33
Equipment.....	33
Chapter 4: Patient Profile	35
Overview	35
How to access Patient Profile.....	35
Summary	35
Flowsheet Calendar – Last 6 Weeks.....	36
Total Alarms & Cautions – Last 20 Treatments	37
Prescription Summary.....	38
Alert Criteria.....	38
Lifestyle Goals	38
Treatments	39
Sort Treatment list.....	39
Therapy Management.....	41
Prescription	41
Therapy Status	43
Alert Criteria.....	44
Chapter 5: Treatment Flowsheet	46
Overview	46
Treatment Status.....	47
How to access the treatment flowsheet	48
Summary	48
Dialysate.....	50
Ultrafiltration.....	50
Blood	51

TOC

Alert, Alarms & Cautions.....	51
Weight.....	52
Prescription.....	53
Clinic Note.....	53
Vitals.....	54
Medication.....	54
Patient Notes & Photos.....	54
Cycler Data.....	56
Vitals, Medication & Equipment.....	58
Vitals.....	58
Medication.....	58
Equipment Data.....	59
Chapter 6: User Management.....	60
Overview.....	60
User Information.....	61
User Search.....	62
Create User.....	63
Description of User Roles.....	64
Edit User Information.....	65



Chapter 1: Indications for Use

Overview

The kinexus Therapy Management Service for HHD is a cloud-based data repository and communication system used between a healthcare facility and chronic home haemodialysis (HHD) patients. The system securely transmits information from patients' NxStage Cyclers via the Nx2me App to the cloud where the data can be stored, aggregated, and presented to healthcare providers in kinexus.

The system is not a medical device and does not control NxStage cyclers. This system shall not be used as a substitute for direct medical intervention or emergency care. Information collected and transmitted requires clinical judgment by an experienced medical professional.

Safety Notice

Note the following:

- The kinexus HHD provider portal does not provide medical diagnoses or recommendations regarding medical treatment.
- The kinexus HHD provider portal and the content appearing on this site are not substitutes for professional medical advice, diagnosis, or treatment.
- In an emergency, patients should call their local emergency assistance number.

System Requirements

Note the following system requirements to run the kinexus HHD provider portal.

Screen Size

The kinexus HHD provider portal is optimised for viewing on laptops and desktop devices with a minimum resolution of 1280x1024.

Operating System and Browser

The kinexus HHD Provider Portal is supported on the latest versions of Microsoft Edge and Google Chrome web browsers.

Setup

Note the following setup information.

Initial Onboarding and Customer Service

The kinexus HHD portal is a web-based application that is typically administered by your organisation. The one-time onboarding process includes setting up your organisation within the kinexus HHD portal and assigning at least one admin user role that will facilitate the roll-out of the portal to your broader home therapy staff and clinics.

Customisation of the Portal

During initial onboarding, some changes may be made to the default settings of the user interface to better suit the needs of an organisation. As a result, some functions described, and screen captures shown, in this IFU may differ slightly from the customised version of the portal that is used at your organisation. See [Support Resources](#) on the next page for more information.

Support Resources

Note the following support resources.

Customer Service

For help with patient enrolment, equipment, and supplies ordering, please contact your regional Fresenius Customer Service office.

Technical Support

For help with troubleshooting or assistance with the kinexus HHD Portal, please contact +44 800 048 8352.

Warnings

A warning alerts the user to the possibility of serious injury, other adverse reactions, or death associated with the use or misuse of the device.

Making clinical decisions



WARNING:

Never use the kinexus HHD Portal as the only method of making clinical decisions, monitoring, or communicating with the patient. Risk of poor therapy outcomes.

The kinexus HHD Portal must be used in conjunction with other methods of assessing clinical signs and symptoms to ensure proper management of patient care.

Manually adjusting cyclers settings



WARNING:

Changes to settings in the kinexus HHD Portal do not update settings on the cycler. The patient or care partner must manually adjust settings on the cycler. Risk of poor therapy outcomes.

Symbols and conventions

This user manual uses a standard set of symbols and formatting conventions.

Symbols and definitions

	WARNING
	NOTE
	Manufacturer

User manual formatting conventions


User interface elements	This formatting convention is used to identify user interface elements such as fields, buttons, tabs, options, and dialog boxes. For example, click OK in the User Assistance dialog box.
--------------------------------	---

Accessing the User Manual

The user manual can be downloaded as a PDF within kinexus HHD.

1. Clicking the question mark (?) located in the navigation pane, [Figure 1-1](#) (A), opens a pop-up window displaying a list of available user manuals, [Figure 1-2](#).


Figure 1-1 Access User Manuals.

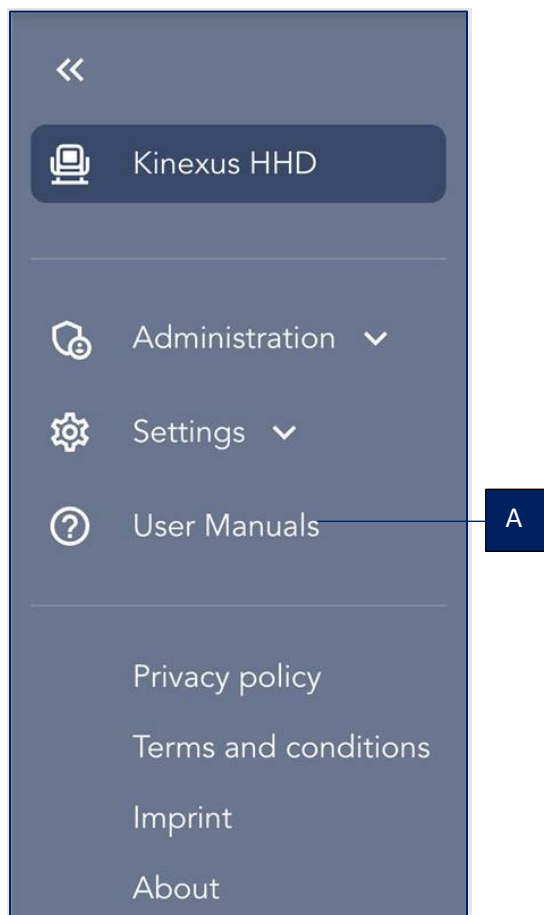

NOTE:

Non-Admin users will not have access to the Administration section in the side navigation menu.

Figure 1-2 List of available User Manuals.

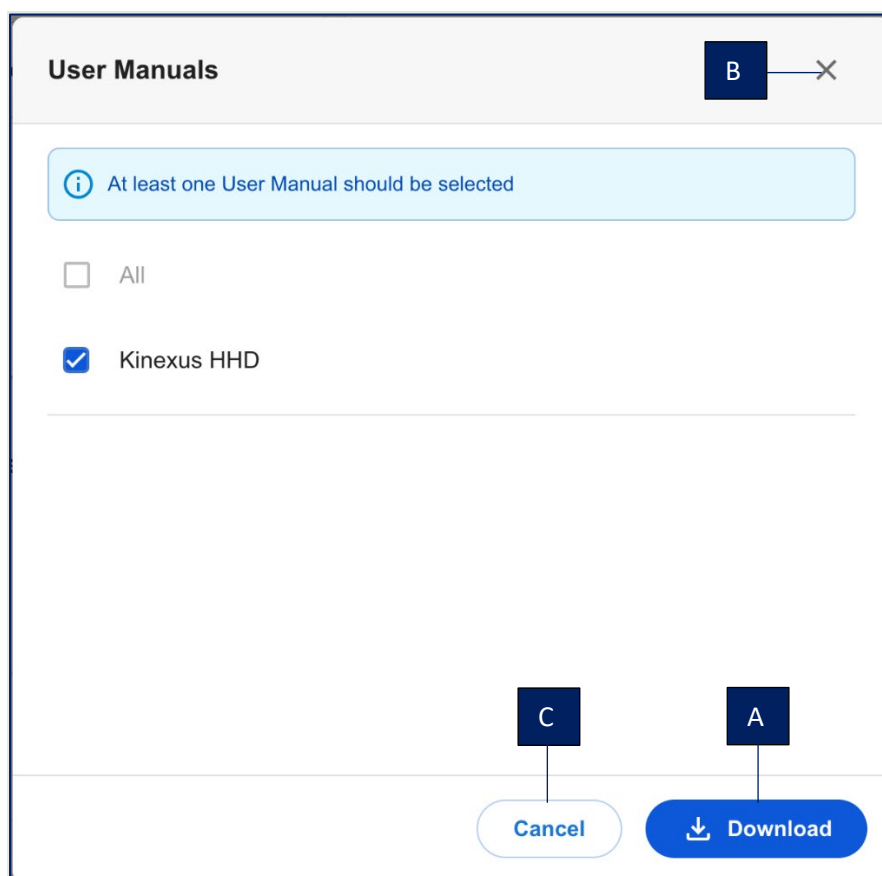


2. Within the User Manuals window, select the check box corresponding to the user manual you would like to download.
3. Click the **Download** button, Figure 1-2 (A), to save the user manual locally for future reference.
4. To close the User Manual window, either click the **X**, Figure 1-2 (B) or select **Cancel** Figure 1-2 (C).

Logging into kinexus HHD Portal

Users will receive the URL from their organisation for logging into the portal. Clicking the URL will redirect the users to the page, as shown in, [Figure 1-3](#). User accounts can be created by admin users only, or users can initiate an account request by selecting **Request to register**, [Figure 1-3](#) (A), and filling out the required information, but access will not be available until an admin approves the request. A valid username and password are required to log into the portal.

Figure 1-3 Login screen.

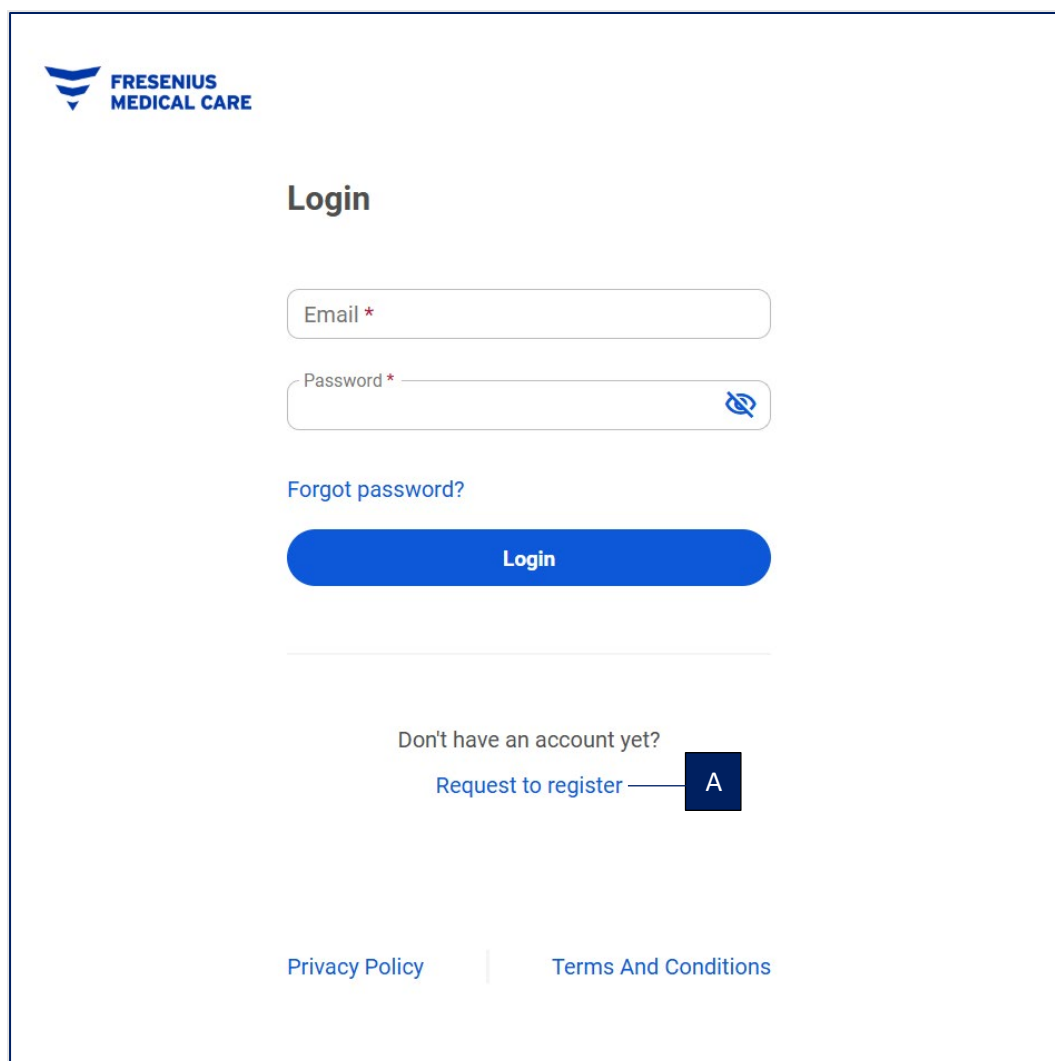

The screenshot shows the login interface for the kinexus HHD Portal. At the top left is the Fresenius Medical Care logo. The main heading is "Login". Below it are two input fields: "Email *" and "Password *". The password field has an eye icon to its right for toggling visibility. Below the password field is a link "Forgot password?". A large blue "Login" button is centered below these fields. A horizontal line separates the login section from the registration section. Below the line, the text "Don't have an account yet?" is followed by a link "Request to register" and a dark blue square button labeled "A". At the bottom, there are two links: "Privacy Policy" and "Terms And Conditions".



NOTE:

After login the user might be required to set up Multi-Factor Authentication (MFA) based on your organisation settings. For more information on MFA please refer to Multi-Factor Authentication (MFA).



NOTE:

Users are required to enter their registered email address; the password field is case-sensitive and masked for security. (It can be unmasked by clicking on the eye symbol next to the field.) Ensure that capitalisation and/or special characters are entered correctly.

Multi-Factor Authentication (MFA)

Multi-Factor Authentication (MFA) is a security process that requires users to provide two or more distinct authentication factors to verify their identity before accessing the kinexus application. MFA involves the generation of a one-time code sent via SMS to the user's mobile number or via email address to the user's email address. The user can then use this code to complete the login process and securely access the application.



NOTE:

MFA Renewal: The MFA mechanism will require renewal under circumstances where the user does not engage in any login activity for a period exceeding 90 days (about 3 months) or opts to modify their web browser configuration. If any of these 2 conditions are not met, MFA will need to be enrolled during their initial login.

Upon user login to the kinexus application, [Figure 1-3](#), a prompt will be displayed to provide your phone number, as shown in [Figure 1-4](#) below.

Figure 1-4 Provide your phone number screen.

Select the appropriate prefix (country-code) from the drop-down list, input a valid phone number, click the **Continue** button, [Figure 1-4](#) (A), to send a one-time code to the phone number provided. Upon receiving the one-time code via SMS, enter the code to successfully log in to the application.

If you prefer to receive the one-time code via email, click the **code to your email** button, [Figure 1-4](#) (B). Once the one-time code is received to the registered email address, enter the code in the specified field and click **Send**, [Figure 1-5](#) (A). This process enables you to successfully log in to the application.

Figure 1-5 Enter code screen.

FRESENIUS MEDICAL CARE

Enter code

Code has been sent to email

Email *
h*****r@yopmail.com

Code *
814593

Should be 6 digits and contain only numbers

Send **A**

Didn't receive a code?
You can request a new code every 1 minute.

Resend to [email](#) or [mobile](#)

B **C**

If the entered code does not match, an error message will be displayed. Users are allowed 5 attempts to enter the correct code. After reaching the maximum number of attempts, the user will be required to wait for a defined period (240 seconds) before being allowed to resend the code.

If the verification code is not received, you can request to send it again to your email address or phone number by clicking **Resend to email**, Figure 1-5 (B), and **Resend to mobile**, Figure 1-5 (C).

Password Reset

From the top navigation menu, clicking on **Change password**, Figure 1-6 (A), will take user to the Change Password screen, as shown in Figure 1-7. Users can reset their account password from this screen.

Figure 1-6 Password reset.

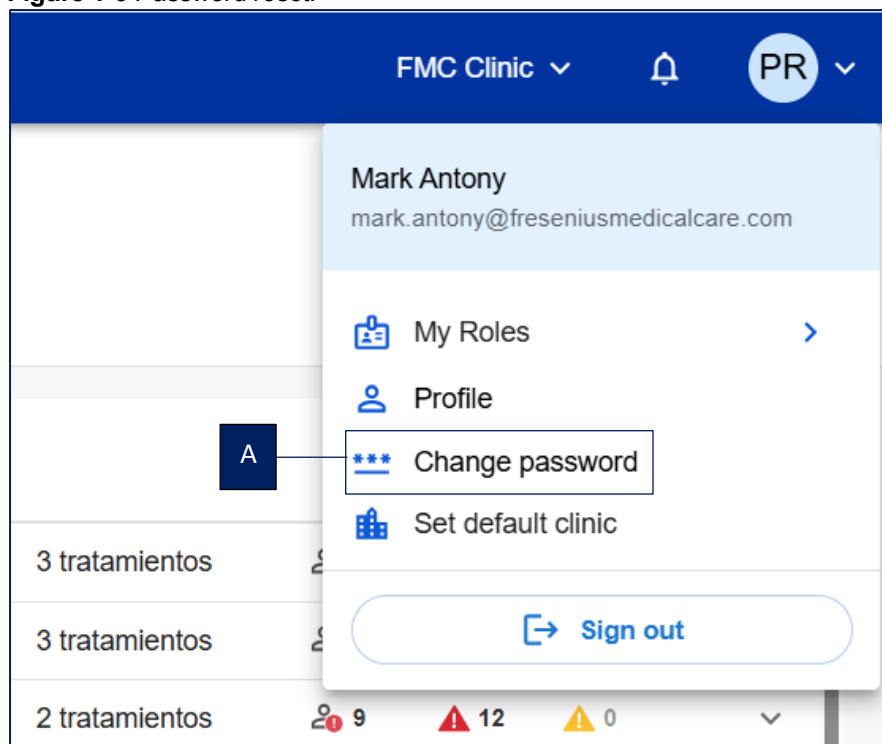


Figure 1-7 Password reset.

FRESENIUS MEDICAL CARE

[< Back](#)

Change password

Current password *

New password *

The password should contain:

- ✓ At least 10 characters
- ✓ A lowercase letter
- ✓ An uppercase letter
- ✓ A number
- ✓ One of the symbols:
! ? \ / , ; < > % # & () * + - = [] { | } ~ ^

Confirm password *

Set Password > **A**

Steps to Reset Your Password

- In the **Current Password** field, enter your existing password.
- In the **New Password** field, type your new password. Ensure that it meets all the requirements listed above, [Figure 1-7](#).
- Re-enter the new password in the **Confirm Password** field to verify it matches the one entered in the New Password field.
- Click the **Set Password** button, [Figure 1-7 \(A\)](#), to save your new password.

If all fields are correctly filled and the new password meets the specified criteria, your password will be successfully updated. If there are any issues, ensure that the new password meets all requirements and that both the New Password and Confirm Password fields match exactly.

END OF CHAPTER

Chapter 2: Introduction to kinexus HHD


Overview

The kinexus HHD Remote Therapy Management Portal helps health care providers monitor patients on home haemodialysis. Patients use electronic rather than paper flowsheets to transfer treatment information to the dialysis centre, as shown in [Figure 2-1](#).

kinexus Therapy Management Service for HHD consists of two main components:

- kinexus HHD Patient Application (Nx2me):: an iOS tablet application (on the app store as Nx2me).
- kinexus HHD Remote Therapy Management Portal (“portal”) : a web-based application for health care providers.

Figure 2-1 *kinexus Therapy Management Service for HHD.*

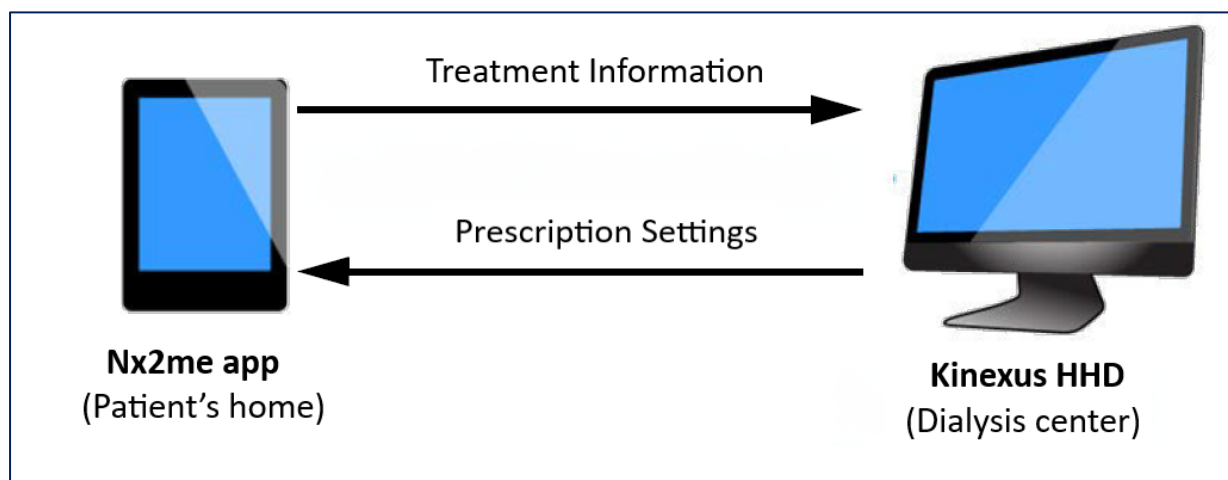

kinexus HHD Patient App (Nx2me)

The patient app is used at home by patients and care partners. The patient app combines treatment information from the cyclor with data entered by the patient or care partner to create a flowsheet that is transferred to the kinexus HHD provider Portal. For more information about patient app, refer to the *kinexus HHD Patient App User Guide (Nx2me)*.

kinexus HHD Provider Portal

The Provider Portal is a web-based application that allows clinicians at the dialysis centre to securely access patient treatment information uploaded from the Nx2me app. Prescription settings for home haemodialysis therapy that are displayed within the Nx2me app are configured within the kinexus HHD Provider Portal

Use

Physicians and clinical staff at dialysis centres use kinexus HHD for remote therapy monitoring to enable the ability to evaluate therapy progress and identify patients that may require intervention.

Disclaimer

The data represented in the screenshots is not medically relevant.

User Roles

The following table provides an overview of the roles available for kinexus HHD users.

User Roles	Rights	Restrictions
User Administrator (Admin)	User can: • View Patient Data. • Create Prescription. • Define Exception Criteria. • View Treatments/Flowsheet Details. • Review/Un-Review Treatments. • Access User Management section.	No restrictions.
Kinexus Home HD User	User can: • View Patient Data. • Create Prescription. • Define Exception Criteria. • View Treatments/Flowsheet Details. • Review/Un-Review Treatments.	User cannot access the User Management section.
Read only User	User can: • View Patient Data. • View Prescription. • View Treatments/Flowsheet Details.	User cannot access the User Management section. User can only view but cannot make any edits.
Patient Administrator	User can: • Add patients to kinexus HHD	This role cannot be used on it's own and must be used in combination with another kinexus HHD access based role. This role may not be available in all regions

END OF CHAPTER


Chapter 3: Clinic Dashboard

Overview

This chapter describes the functions that are available within the kinexus HHD clinic dashboard. The Clinic Dashboard is the landing page upon logging into the portal.

The Clinic Dashboard includes the following tabs to easily access its main functions, [Figure 3-1 \(A\)](#).


Figure 3-1 Navigation tabs in Clinic Dashboard.

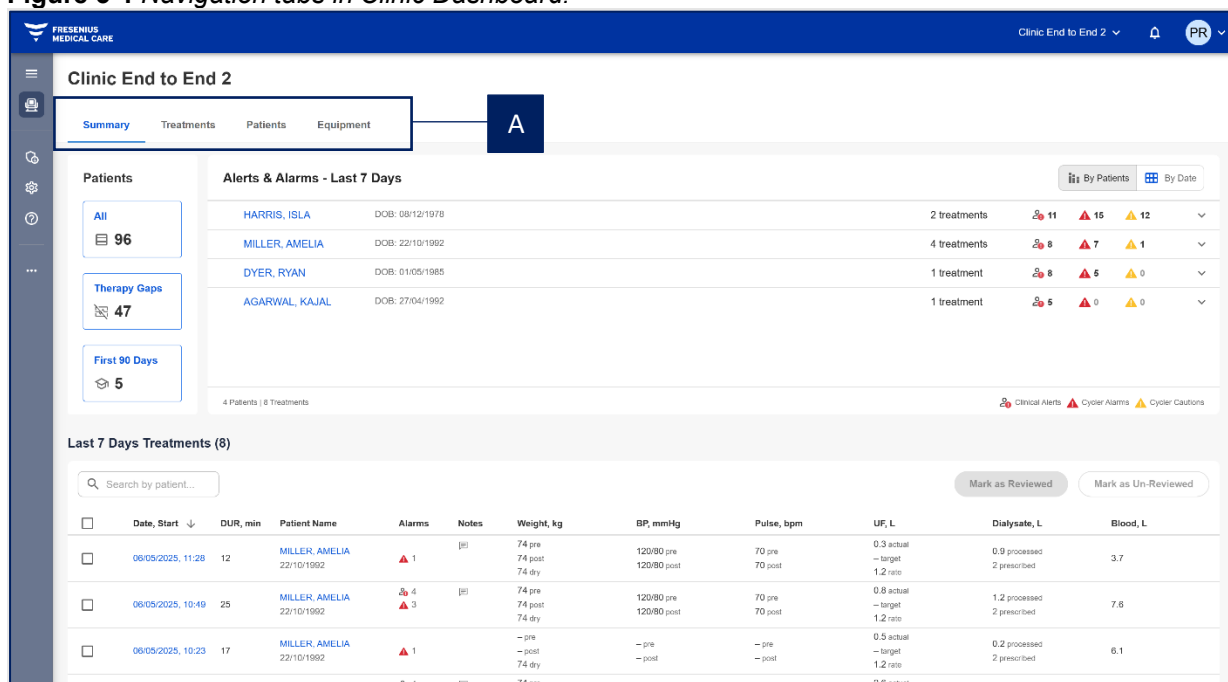

- **Summary** – see [Summary](#) for more information.
- **Treatments** – see [Treatments](#) for more information.
- **Patients** – see [Patients](#) for more information.
- **Equipment** – see [Equipment](#) for more information.

Summary

Figure 3-2 Clinic Dashboard.

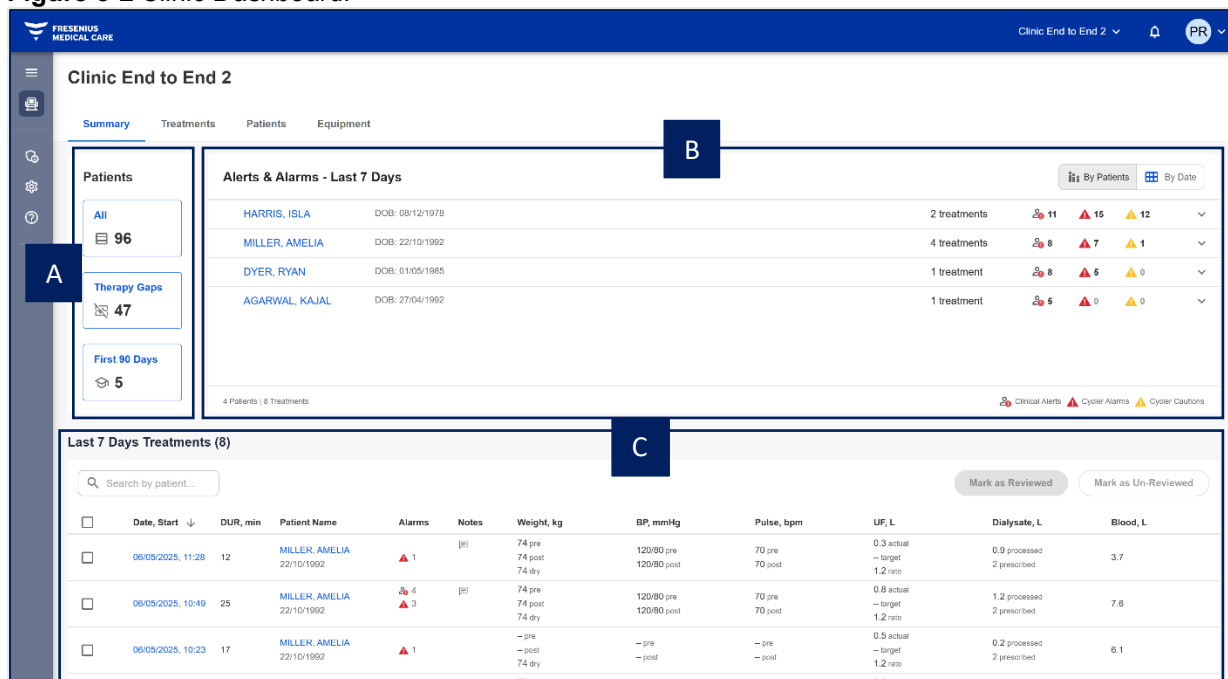


The *Clinic Dashboard – Summary* tab contains three sections:

- **Clinic Overview**, Figure 3-2 (A), provides an overview of **All** patients associated with the clinic and displays the number of patients that currently have **Therapy Gaps**. And the number of patients in their **First 90 Days**, calculated from the date of their first treatment recorded in the kinexus HHD system.
- **Alerts & Alarms – Last 7 Days**, Figure 3-2 (B), displays the patients list along with treatments within the last 7 days that have received **Clinical Alerts**, **Cycler Alarms** and/or **Cycler Cautions**.
- **All Treatments**, Figure 3-2 (C), this table lists flowsheets that the clinic has received for the patients within the clinic for the last 7 days treatments.

Access *Clinic Overview Card*

The *Clinic Overview Card* displays **All Patients**, **Therapy Gaps** and **First 90 days** with a count of the total number of clinic patients associated with each classification. The **All Patients** count includes all clinic patients that are enrolled in kinexus HHD. The **Therapy Gap** count includes all patients that have not sent a treatment within the last 3 days or more. **First 90 Days** count includes all clinic patients calculated from the date of their first treatment recorded in the kinexus HHD system.

Figure 3-3 *Clinic Overview Card.*

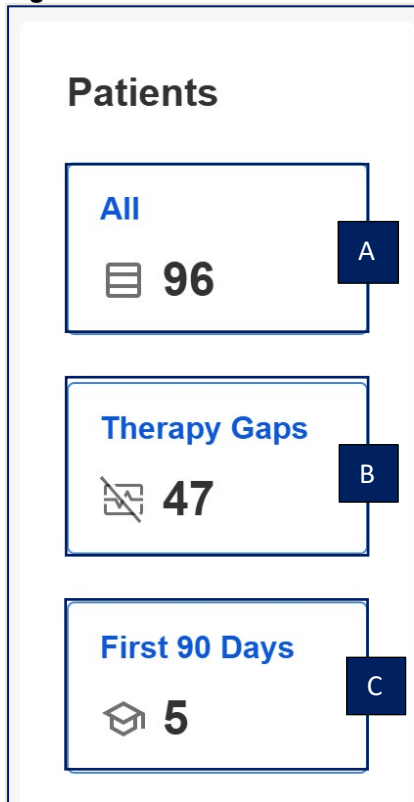

All

Clicking **All** [Figure 3-3 \(A\)](#), takes you to the **Patients** section which displays a list of all active patients in the selected clinic as shown in [Figure 3-4](#). Patients are sorted alphabetically by last name and include details such as HHD status, ID, Sex, date of birth (DOB), Therapy Gap, HHD Start date, 1st treatment date, and the number of unreviewed treatments.

Figure 3-4 Clinic Overview Card.

Summary Treatments (+ 149) Patients Equipment								
Search by patient...		All (293)		+ Add Patient				
Name ↑	HHD Status	ID	Sex	DOB	Therapy Gap	HHD Start	1st Treatment	Un-Reviewed
Archer, Marie	Enrolled	4597563287	Male	01/02/1985		08/05/2025		
Antony, Mark	Enrolled		Male	05/05/1968		21/04/2025		
Archer, Annie	Enrolled		Male	05/05/1987		21/04/2025		
Arden, Ashley	Enrolled	7899965164	Female	05/09/1990		21/04/2025		
Avery, Addison	Enrolled	7733589628	Other	05/05/1988		21/04/2025		
Andrea, Araminta	Enrolled		Male	23/02/1989		23/07/2025		
Alfred, Ashton	Enrolled	25897469625	Male	02/02/1985		29/04/2025		
Augustine, Arnold	Enrolled		Female	20/09/1986		14/08/2025		
Abel, Tesfaye	Enrolled	7853492790	Male	16/05/1983	72 Days	16/06/2025	16/06/2025	
ADAMS, JACKSON	Enrolled		Male	07/03/1988	15 Days	03/07/2024	15/08/2024	3 treatments

Therapy Gaps

Clicking **Therapy Gaps** [Figure 3-3 \(B\)](#), takes you to the **Patients** list, displaying patients who have not recorded a treatment in the past three days or more. For more information, see the [Patients](#) section.

First 90 Days

Clicking **First 90 Days** [Figure 3-3 \(B\)](#), takes you to the **Patients** list, displaying only those patients whose first treatment has been recorded within the last 90 days in the kinexus HHD system. For more information, see the [Patients](#) section.

Alerts & Alarms – Last 7 Days

The Clinical Alerts & Alarms card displays treatments that received a clinical alert and/or cyclus alarm and cautions within the last 7 days. Each row contains the patient's name, DOB, the number of treatments impacted, the count of clinical alerts, the total count of cyclus alarms and cautions. For more information, see the [Alert Criteria](#) section in [Chapter 4: Patient Profile](#).

Clicking the arrow adjacent to the patient's name, [Figure 3-5 \(A\)](#), displays list of treatment dates with the total treatment duration and number of clinical alerts and/or cyclus alarms and cautions, as shown in [Figure 3-6 \(A\)](#).

Figure 3-5 Clinical Alerts and Cyclus Alarms card.

Alerts & Alarms - Last 7 Days					
		By Patients		By Date	
• AJOHN, ADAM	DOB: 15/05/1985	8 treatments	14	89	7
YOUNG, POPPY	DOB: 06/07/1987	4 treatments	10	9	8
• CLARK, JACOB	DOB: 24/04/1994	2 treatments	10	11	
David, Tim	DOB: 30/05/1996	1 treatment	7	2	
• GEMMA, DANIEL	DOB: 28/05/1984	2 treatments	4		
• LAURA, COOPER	DOB: 22/04/1999	2 treatments	2		
• DENNIS, STEVE	DOB: 29/04/1976	4 treatments	0	4	
7 Patients 20 Treatments					
Clinical Alerts Cyclus Alarms Cyclus Cautions					

Figure 3-6 Clinical Alerts and Cyclus Alarms and Caution list.

Alerts & Alarms - Last 7 Days					
		By Patients		By Date	
• AJOHN, ADAM	DOB: 15/05/1985	8 treatments	14	89	7
04/06/2025	07:32, 0 mins	1			
03/06/2025	13:12, 15 mins	2	5	4	
03/06/2025	09:19, 29 mins	2	74	1	
03/06/2025	08:28, 4 mins	3	5		
03/06/2025	07:16, 15 mins	2	4	2	
30/05/2025	10:49, 0 mins		1		
7 Patients 20 Treatments					
Clinical Alerts Cyclus Alarms Cyclus Cautions					

Clicking on the patient's name will take you to the Patient Profile. For more information, see [Chapter 4: Patient Profile](#).

Clicking on the treatments will take you to that treatment flowsheet screen. For more information, see [Chapter 5: Treatment Flowsheet](#).

Clicking **By Patients**, [Figure 3-6 \(B\)](#) displays the alerts and alarms grouped by patient in the table. Clicking **By Date**, [Figure 3-6 \(C\)](#) lists the most recent alerts and alarms received, sorted by date.

Last 7 Days Treatments

The **Last 7 Days Treatments** table lists all treatments recorded in the past seven days for patients associated with the clinic, as shown in [Figure 3-7](#). Treatments are listed chronologically by the newest date. To sort treatments from oldest to newest dates, click the downward arrow next to the **Date, Start** column header, [Figure 3-7](#) (A).

Search by patient

The search bar at the left side of this screen enables users to search for patients by first name or last name. Searching happens automatically once you enter the first letter into the search bar, [Figure 3-7](#) (B). To return to the full patient list, remove your search parameters from the search bar.

Figure 3-7 Last 7 Days Treatments table.

Last 7 Days Treatments (34)											
Search by patient...											
Selected: 2											
☐	Date, Start	DUR, min	Patient Name	Alarms	Notes	Weight, kg	BP, mmHg	Pulse, bpm	UF, L	Dialysate, L	Blood, L
☐	05/06/2025, 06:41	0	LEWIS, EMILY 08/10/1982	4		40 pre 40 post 40 dry	120/80 pre 120/80 post	70 pre 70 post	-- actual 2.2 target -- rate	0 processed 6 prescribed	
☒	05/06/2025, 06:20	0	LAURA, COOPER 22/04/1999	1		-- pre -- post 60.5 dry	-- pre -- post	-- pre -- post	-- actual -- target -- rate	0 processed 5.6 prescribed	
☒	05/06/2025, 07:55	32	LEWIS, EMILY 08/10/1982	5 1 2		45.1 pre 45.1 post 40 dry	110/80 pre 110/110 post	70 pre 85 post	1 actual 4.2 target 1.6 rate	3 processed 6 prescribed	6.7
☐	05/06/2025, 07:49	0	LAURA, COOPER 22/04/1999	2		-- pre -- post 60.5 dry	-- pre -- post	-- pre -- post	-- actual -- target -- rate	0 processed 5.6 prescribed	
☐	04/06/2025, 12:23	0	YOUNG, POPPY 06/07/1987	12		62.8 pre 62.8 post 62.8 dry	120/116 pre -- post	70 pre -- post	-- actual -- target -- rate	0 processed 20 prescribed	
☐	04/06/2025, 11:43	12	YOUNG, POPPY 06/07/1987	6 1		62.8 pre 62.8 post 62.8 dry	139/93 pre 139/93 post	70 pre 70 post	0.2 actual 1 target 0.43 rate	0.2 processed 20 prescribed	2.8
☐	04/06/2025, 10:39	15	YOUNG, POPPY 06/07/1987	1 1		64.8 pre 62.8 post 62.8 dry	120/80 pre 120/80 post	70 pre 70 post	0.3 actual 2 target -- rate	0.6 processed 20 prescribed	3.5

Last 7 Days Treatments Table

The following information is displayed in the *Last 7 Days Treatments* table, [Figure 3-7](#) (C):

- **Date, Start** displays the date when the treatment started in DD/MM/YY format. A blue dot beside the date means the treatment has not yet been reviewed. Clicking on the treatment date will take you to that patient's treatment flowsheet page. For more information, see [Chapter 5: Treatment Flowsheet](#).
- **DUR, min** displays the total length of the treatment in minutes.
- **Patient Name** displays the patient's name (last name, first name) and date of birth. Clicking on the patient's name will take you to that patient profile page. For more information, see [Chapter 4: Patient Profile](#).
- **Alarms** indicators appear when a treatment value falls outside a prescribed range as defined by the alert criteria.

 1 The number adjacent to the alert icon indicates the total count of clinical alerts.

 12 The number adjacent to the alert icon indicates the total count of cyclor alarms.

 1 The number adjacent to the alert icon indicates the total count of cyclor cautions.

- **Note** indicators appear where the patient/clinic attaches their photos or notes related to the treatment.



A message icon indicates that a text note is attached.



A camera icon indicates that a photo is attached.



A clinic note icon indicates that a text note is attached from the clinic.

- **Weight, kg** displays the patient's pre and post dialysis weight along with the dry weight in kilogrammes (kg). Dry weight is defined in the patient's prescription therapy.
- **BP, mmHg** displays the patient's pre and post blood pressure recorded during the treatment.
- **Pulse, bpm** displays the patient's pre and post pulse recorded during the treatment.
- **UF, L** displays the actual and target Ultrafiltration volume of the treatment along with rate, in litres (L). Rate is the cyclor set rate of ultrafiltration.
- **Dialysate, L** displays the prescribed and processed amount of dialysate used during the entire treatment, in litres (L).
- **Blood, L** displays the total amount of blood processed during the entire treatment, in litres (L).

How to Review Treatments

Clicking the checkbox next to a specific treatment date, [Figure 3-8 \(A\)](#), enables the **Mark as Reviewed** or **Mark as Un-Reviewed** buttons, [Figure 3-8 \(B\)](#), at the top of the table.

Users should follow clinic policy and procedure regarding what constitutes a “Review” when utilising this feature.

Figure 3-8 Review Treatments.

Last 7 Days Treatments (34)												
Search by patient... B Mark as Reviewed Mark as Un-Reviewed												
Selected: 2												
☐	Date, Start	DUR, min	Patient Name	Alarms	Notes	Weight, kg	BP, mmHg	Pulse, bpm	UF, L	Dialysate, L	Blood, L	
☐	05/06/2025, 08:41	0	LEWIS, EMILY 08/10/1982	4		40 pre 40 post 40 dry	120/80 pre 120/80 post	70 pre 70 post	— actual 2.2 target — rate	0 processed 6 prescribed		
☒	05/06/2025, 08:20	0	LAURA, COOPER 22/04/1999	1		— pre — post 60.5 dry	— pre — post	— pre — post	— actual — target — rate	0 processed 5.6 prescribed		
☒	05/06/2025, 07:55	32	LEWIS, EMILY 08/10/1982	5 1 2		45.1 pre 45.1 post 40 dry	110/80 pre 110/110 post	70 pre 85 post	1 actual 4.2 target 1.6 rate	3 processed 6 prescribed	6.7	
☐	05/06/2025, 07:49	0	LAURA, COOPER 22/04/1999	2		— pre — post 60.5 dry	— pre — post	— pre — post	— actual — target — rate	0 processed 5.6 prescribed		
☐	04/06/2025, 12:23	0	YOUNG, POPPY 06/07/1987	12		62.8 pre 62.8 post 62.8 dry	120/116 pre — post	70 pre — post	— actual — target — rate	0 processed 20 prescribed		
☐	04/06/2025, 11:43	12	YOUNG, POPPY 06/07/1987	6 1		62.8 pre 62.8 post 62.8 dry	139/93 pre 139/93 post	70 pre 70 post	0.2 actual 1 target 0.43 rate	0.2 processed 20 prescribed	2.8	
☐	04/06/2025, 10:39	15	YOUNG, POPPY 06/07/1987	1 1		64.8 pre 62.8 post 62.8 dry	120/80 pre 120/80 post	70 pre 70 post	0.3 actual 2 target — rate	0.6 processed 20 prescribed	3.5	

The **Mark as Reviewed** button only appears if you select a treatment that has not been reviewed; the **Mark as Un-Reviewed** button only appears if you select a treatment that has been reviewed. You may select multiple treatments at a time by clicking checkboxes, or you may select all displayed treatments by clicking the checkbox next to Date, Start at the top of the table.

If you select both **Reviewed** and **Unreviewed** treatments, both the **Mark as Reviewed** and **Mark as Un-Reviewed** buttons will appear.

Clicking on the patient's name will take you to the Patient Profile. For more information, see [Chapter 4: Patient Profile](#).

Clicking on treatment will take you to the Treatment Flowsheet. For more information, see [Chapter 5: Treatment Flowsheet](#).

Treatments

The **Treatments** tab, [Figure 3-9](#), lists treatments according to start date and time. Only treatments that your user role and designated clinic authorize you to view will be visible. To sort treatments from oldest to newest dates, click the downward arrow next to the **Date, Start** column header. For more information, see section, [Last 7 Days Treatments](#).

Figure 3-9 *Treatments tab.*


Summary Treatments (• 149) Patients Equipment											
A		B		D							
Q Search by patient...		All (2083)		Mark as Reviewed Mark as Un-Reviewed							
☐	Date, Start ↓	DUR, min	Patient Name	Alarms	Notes	Weight, kg	BP, mmHg	Pulse, bpm	UF, L	Dialysate, L	Blood, L
☐	10/10/2025, 17:19	33	GEMMA, DANIEL 28/05/1984	2 4 3	— pre — post 39.9 dry	— pre — post	— pre — post	— pre — post	1.1 actual — target 1.57 rate	1.9 processed 56 prescribed	10.9
☐	10/10/2025, 13:37	43	GEMMA, DANIEL 28/05/1984	2 6 4	— pre — post 39.9 dry	— pre — post	— pre — post	— pre — post	1.4 actual — target 2.21 rate	1.5 processed 56 prescribed	14.3
☐	09/10/2025, 21:08	0	Ciphra, Tony 19/02/1999	1	— pre — post 70 dry	— pre — post	— pre — post	— pre — post	— actual — target — rate	0 processed 2.3 prescribed	
☐	09/10/2025, 17:11	0	Ciphra, Tony 19/02/1999	1	— pre — post 70 dry	— pre — post	— pre — post	— pre — post	— actual — target — rate	0 processed 2.3 prescribed	
☐	09/10/2025, 08:33	0	GEMMA, DANIEL 28/05/1984		— pre — post 39.9 dry	— pre — post	— pre — post	— pre — post	— actual — target — rate	0 processed 56 prescribed	
☐	09/10/2025, 01:07	0	Ciphra, Tony 19/02/1999		— pre — post 70 dry	— pre — post	— pre — post	— pre — post	— actual — target — rate	0 processed 2.3 prescribed	
☐	09/10/2025, 00:36	0	Ciphra, Tony 19/02/1999	1	— pre — post 70 dry	— pre — post	— pre — post	— pre — post	— actual — target — rate	0 processed 2.3 prescribed	

The **Search by patient**, [Figure 3-9](#) (A), bar at the top left side of this screen enables you to search for patients by first name or last name. For more information, refer to the [Search by patient](#) section.

Sort Treatment list

At the top of the *Treatments* list is a sort drop-down menu, [Figure 3-9](#) (B), which allows you to sort the list in four ways, as shown in [Figure 3-10](#).

Figure 3-10 Sort treatment drop-down list.

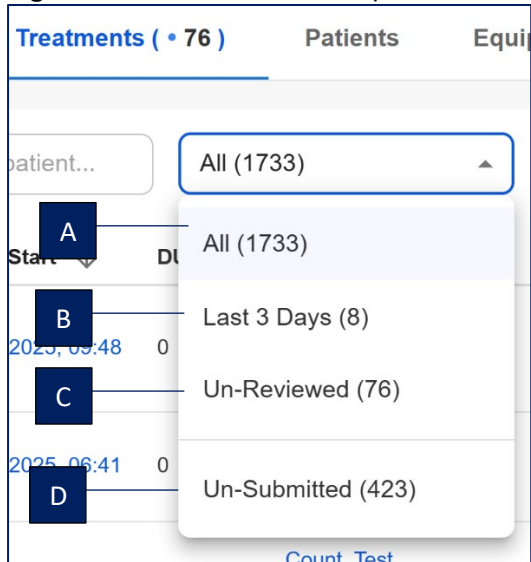

- **All**, [Figure 3-10](#) (A), displays a complete list of treatments received by all patients associated with the clinic.
- **Last 3 days**, [Figure 3-10](#) (B), displays treatments received within the last three days.
- **Un-Reviewed**, [Figure 3-10](#) (C), displays treatments that have not yet been reviewed.
- **Un-Submitted**, [Figure 3-10](#) (D), displays treatments that patients have not yet submitted and remain editable.

The treatment table, [Figure 3-9](#) (C), displays a list of information related to each treatment. For more information on the treatment table, see the [Last 7 Days Treatment Table](#) section.

Clicking the checkbox next to a specific treatment date, enables the **Mark as Reviewed** or **Mark as Un-Reviewed**, [Figure 3-9](#) (D), buttons, at the top of the table. For more information. See [How to Review Treatments](#) section.

Clicking on the patient's name will take you to the Patient Profile screen. For more information, see [Chapter 4: Patient Profile](#).

Clicking on treatment will take you to the Treatment Flowsheet screen. For more information, see [Chapter 5: Treatment Flowsheet](#).

Patients

The **Patients** tab, [Figure 3-11](#), displaying a list of all active patients associated with that specific clinic. Patients are listed alphabetically by last name, and first name. Additionally, a search field is provided to facilitate searching for specific patients within the list. To sort patients by last name, click the upward arrow next to the **Name** column header.

Figure 3-11 *Patients tab.*


Summary Treatments (+ 149) Patients Equipment									
A Search by patient...		B All (293)					D + Add Patient		
Name ↑	HHD Status	ID	Sex	DOB	Therapy Gap	C	HHD Start	1st Treatment	Un-Reviewed
Archer, Marie	Enrolled	4597563287	Male	01/02/1985			08/05/2025		⋮
Antony, Mark	Enrolled		Male	05/05/1968			21/04/2025		⋮
Archer, Annie	Enrolled		Male	05/05/1987			21/04/2025		⋮
Arden, Ashley	Enrolled	7899965164	Female	05/09/1990			21/04/2025		⋮
Avery, Addison	Enrolled	7733589628	Other	05/05/1988			21/04/2025		⋮
Andrea, Araminta	Enrolled		Male	23/02/1989			23/07/2025		⋮
Alfred, Ashton	Enrolled	25897469825	Male	02/02/1985			29/04/2025		⋮
Augustine, Arnold	Enrolled		Female	20/09/1986			14/08/2025		⋮
Abel, Tesfaye	Enrolled	7853492790	Male	16/05/1983	72 Days		16/06/2025	16/06/2025	⋮
ADAMS, JACKSON	Enrolled		Male	07/03/1988	15 Days		03/07/2024	15/08/2024	3 treatments ⋮

The **Search by patient**, [Figure 3-11](#) (A), bar at the top left side of this screen enables you to search for patients by first name or last name. For more information, refer to the [Search by patient](#) section.

Sort Patients list

At the top of the *Patients list* is a sort drop-down menu, [Figure 3-11](#) (B), which allows you to sort the list in four ways, as shown in [Figure 3-12](#).

Figure 3-12 Sort patients drop-down list.

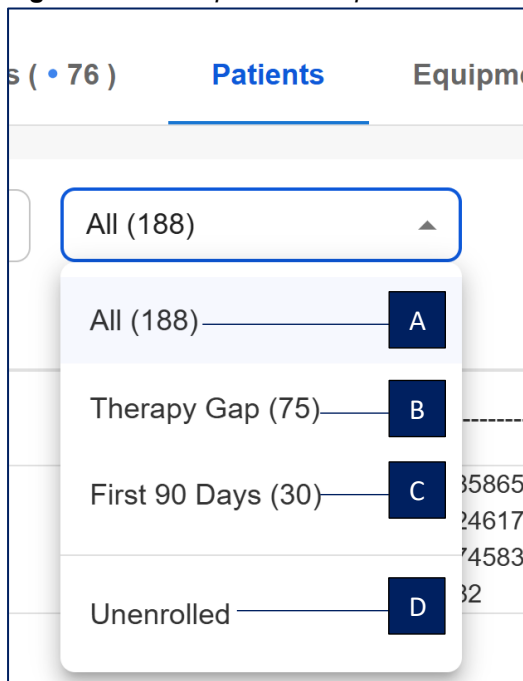

- **All**, [Figure 3-12](#) (A), displays a list of all active patients associated in the selected clinic.
- **Therapy Gaps**, [Figure 3-12](#) (B), displays a list of patients who have not recorded a treatment in the past three days or more.
- **First 90 Days**, [Figure 3-12](#) (C), displays only those patients whose first treatment has been recorded within the last 90 days in the kinexus HHD system.
- **Unenrolled**, [Figure 3-12](#) (D), displays a list of patients who have been unenrolled from the kinexus HHD system.

Patients table

The Patients table displays detailed information about each patient associated with the clinic, including their personal details, clinic information, historical treatment data, and any therapy programs managed within the kinexus HHD system.

The following information is displayed in the *Patients* table, [Figure 3-11](#) (C):

- **Name** is displayed in Last Name, First Name order.
- **HHD Status** displays the patient's current status in the kinexus HHD system. Patients may have statuses such as **Enrolled**, **First 90 Days**, or **Unenrolled**
- **ID** is the identification number that is created by the patient's healthcare provider.

- **Sex** displays the patient's biological designation.
- **DOB** displays the patient's date of birth.
- **Therapy Gap** is only displayed if it has been more than 3 days since the patient's last treatment, it displays the number of days since the patient's last recorded treatment.

HHD Start displays the date when the patient was enrolled in HHD therapy (as entered by the provider).

- **1st Treatment** displays the date of the patient's first treatment recorded in the kinexus HHD system.
- **Un-Reviewed** displays the total number of treatments that has not yet been reviewed in the kinexus HHD system.


+ Add Patient

Clicking the **+ Add Patient** [Figure 3-11](#) (D), button opens a multi-step form that allows you to fill out information to register a patient with your organization. To complete the patient enrollment process, enter the appropriate information in the **Create patient** form.

Equipment

The **Equipment** tab displays detailed information about the equipment used in the patient's treatment, as shown in [Figure 3-13](#).

Figure 3-13 *Equipment tab.*

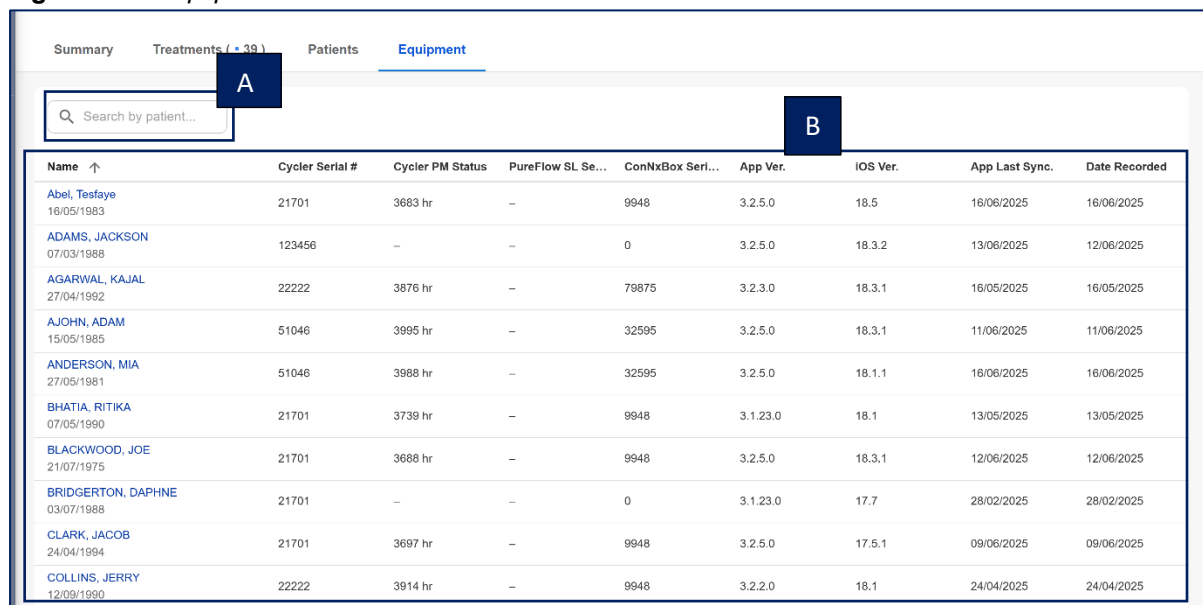

Name ↑	Cycler Serial #	Cycler PM Status	PureFlow SL Se...	ConNxBox Seri...	App Ver.	IOS Ver.	App Last Sync.	Date Recorded
Abel, Tesfaye 16/05/1983	21701	3683 hr	—	9948	3.2.5.0	18.5	16/06/2025	16/06/2025
ADAMS, JACKSON 07/03/1988	123456	—	—	0	3.2.5.0	18.3.2	13/06/2025	12/06/2025
AGARWAL, KAJAL 27/04/1992	22222	3876 hr	—	79875	3.2.3.0	18.3.1	16/05/2025	16/05/2025
AJOHN, ADAM 15/05/1985	51046	3995 hr	—	32595	3.2.5.0	18.3.1	11/06/2025	11/06/2025
ANDERSON, MIA 27/05/1981	51046	3988 hr	—	32595	3.2.5.0	18.1.1	16/06/2025	16/06/2025
BHATIA, RITIKA 07/05/1990	21701	3739 hr	—	9948	3.1.23.0	18.1	13/05/2025	13/05/2025
BLACKWOOD, JOE 21/07/1975	21701	3688 hr	—	9948	3.2.5.0	18.3.1	12/06/2025	12/06/2025
BRIDGERTON, DAPHNE 03/07/1988	21701	—	—	0	3.1.23.0	17.7	28/02/2025	28/02/2025
CLARK, JACOB 24/04/1994	21701	3697 hr	—	9948	3.2.5.0	17.5.1	09/06/2025	09/06/2025
COLLINS, JERRY 12/09/1990	22222	3914 hr	—	9948	3.2.2.0	18.1	24/04/2025	24/04/2025

The **Search by patient**, [Figure 3-13](#) (A), bar at the top left side of this screen enables you to search for patients by first name or last name. For more information, refer to the [Search by patient](#) section.

The following information is displayed in the Equipment table, [Figure 3-13 \(B\)](#):

- **Name** displays the patient's last name followed by their first name, along with their date of birth (DOB) in DD/MM/YY format.
- **Cycler Serial #**, displays the serial number of the cycler as recorded by the last treatment.
- **Cycler PM Status** indicates the remaining cycler usage hours before preventive maintenance is required.
- **PureFlow SL Serial #**, displays the serial number of PureFlow unit as recorded by the last treatment.
- **ConNxBox Serial #**, displays the serial number of the ConNxBox as recorded by the last treatment.
- **App Ver** displays the version of kinexus HHD Patient App (also known as Nx2me).
- **iOS Ver** displays the iOS version of the iPad that is being used with the Patient App.
- **App Last Sync** displays the date when the app was last connected to the internet and synced with the kinexus HHD portal.
- **Date Recorded**, displays the date when the ConNxBox, PureFlow, and Cycler Serial numbers was recorded or gathered.

Clicking on the patient's name will take you to the Patient Profile screen. For more information, see [Chapter 4: Patient Profile](#).

END OF CHAPTER

Chapter 4: Patient Profile

Overview

This chapter provides detailed information about the data available within the kinexus HHD patient profile, including:

- All treatment data
- Prescription details
- Alert criteria
- Lifestyle goals
- Flowsheet calendar report
- Trending alarm & caution data

How to access Patient Profile

Clicking on any of the patient's names, will display the respective patient's profile screen. Here, you can view the patient's demographic details, including date of birth (DOB), Sex, unique identifiers (such as NHS number), and HHD Start Date.

The Patient Profile screen includes the following tabs to easily access its main functions.


- **Summary** – see [Summary](#) for more information.
- **Treatments** – see [Treatments](#) for more information.
- **Therapy Management** – see [Therapy Management](#) for more information.

Summary

The **Patient Profile – Summary** tab contains the following details:

- **Flowsheet Calendar - Last 6 Weeks**, [Figure 4-1](#) (A)
- **Total Alarms & Cautions - Last 20 Treatments**, [Figure 4-1](#) (B)
- **Prescription Summary**, [Figure 4-1](#) (C)
- **Alert Criteria**, [Figure 4-1](#) (D)
- **Lifestyle Goals**, [Figure 4-1](#) (E)

Figure 4-1 Patient profile – Summary tab.

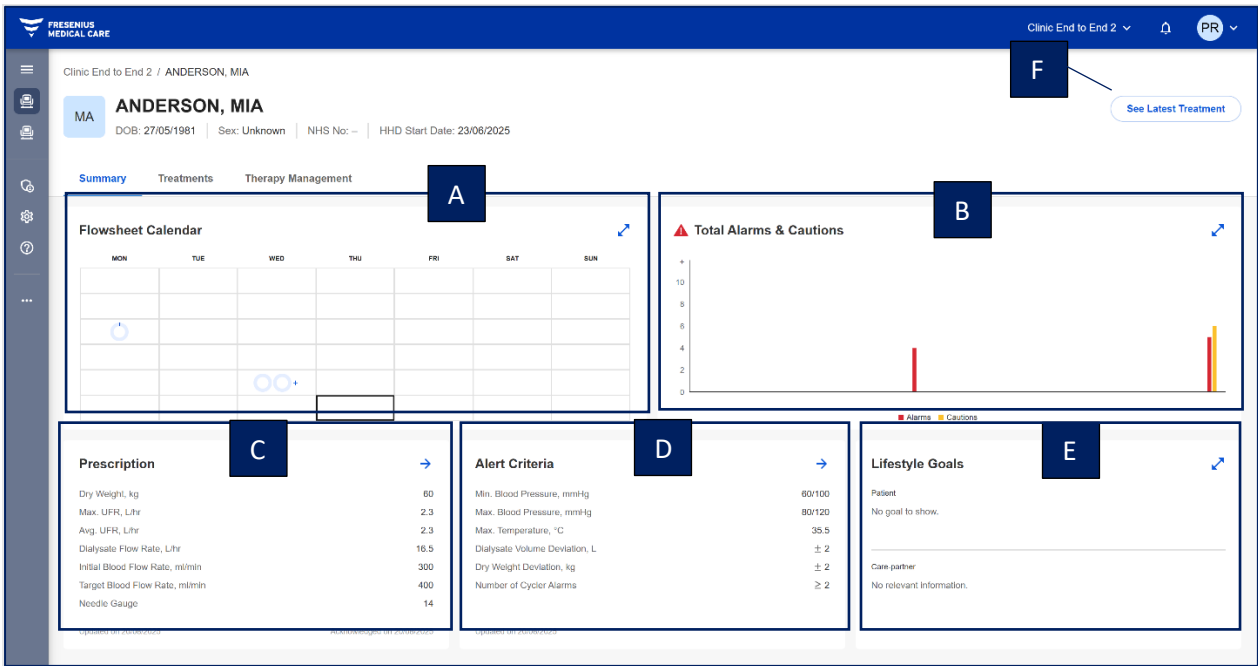


Flowsheet Calendar – Last 6 Weeks

The Flowsheet Calendar displays the percentage of prescribed dialysate that was processed during the treatment (percentage = dialysate processed / dialysate prescribed). The calendar will display this information for all treatments received during the last six calendar weeks, including both confirmed and unconfirmed flowsheets.

Clicking the  icon, Figure 4-1 (A), opens a pop-up window, as shown in Figure 4-2.

Figure 4-2 Flowsheet Calendar.

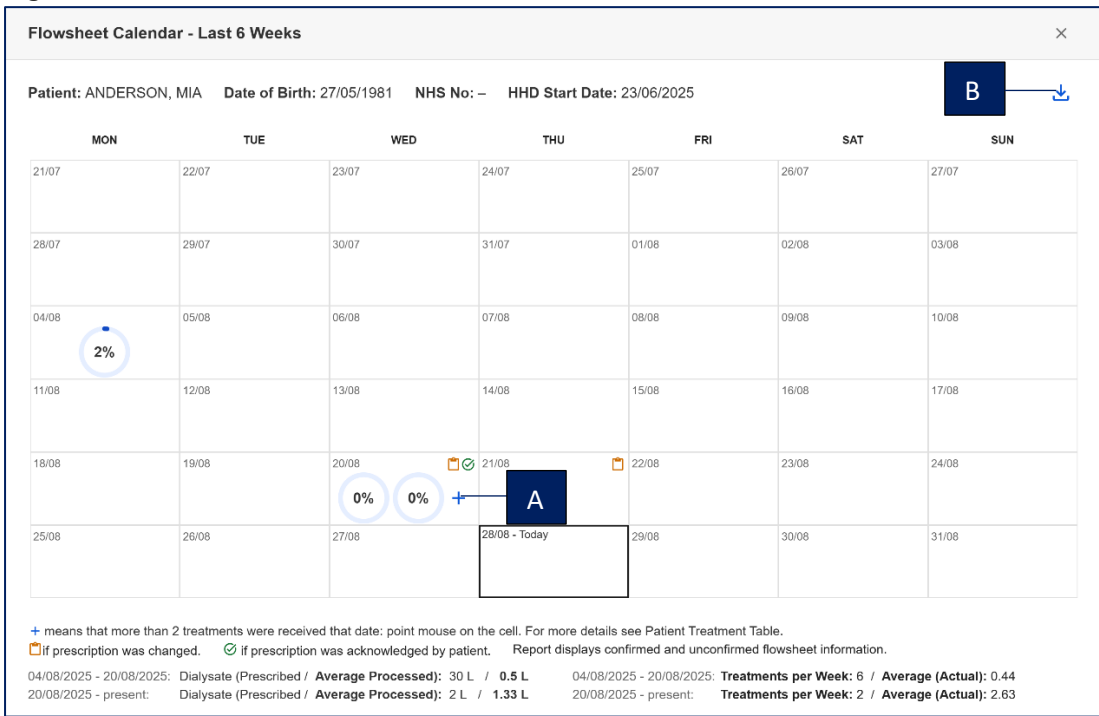

If the patient performed two treatments on the same day, two graphics will appear for that day. If the patient performed more than two treatments in the same day, the + icon, [Figure 4-2 \(A\)](#), appears, in this scenario please refer to the Treatment Table to view information for all treatments.

Hovering your mouse over the graphics will display the treatment date, time, and the percentage of dialysate completed.

Clicking on graphics will take you to the Treatment Flowsheet screen. For more information, see [Chapter 5: Treatment Flowsheet](#).


Clicking on download icon, [Figure 4-2 \(B\)](#), allows you to download and save the Flowsheet Calendar as a PDF, which can then be printed on a local printer.

Total Alarms & Cautions – Last 20 Treatments

Any alarms or cautions that occur on the NxStage cyclor during dialysis are recorded in the treatment flowsheet. The *Total Alarms & Cautions* card provides a cross-treatment view, showing the total number of alarms and cautions across all treatments. The graph displays the number of alarms and/or cautions that occurred during the last 20 treatments. If the number of alarms and/or cautions exceeds ten, they are represented as "10+."

Clicking the [↗](#) icon, [Figure 4-1 \(B\)](#), opens a pop-up window, as shown in [Figure 4-3](#).

Figure 4-3 *Total Alarms & Cautions screen.*

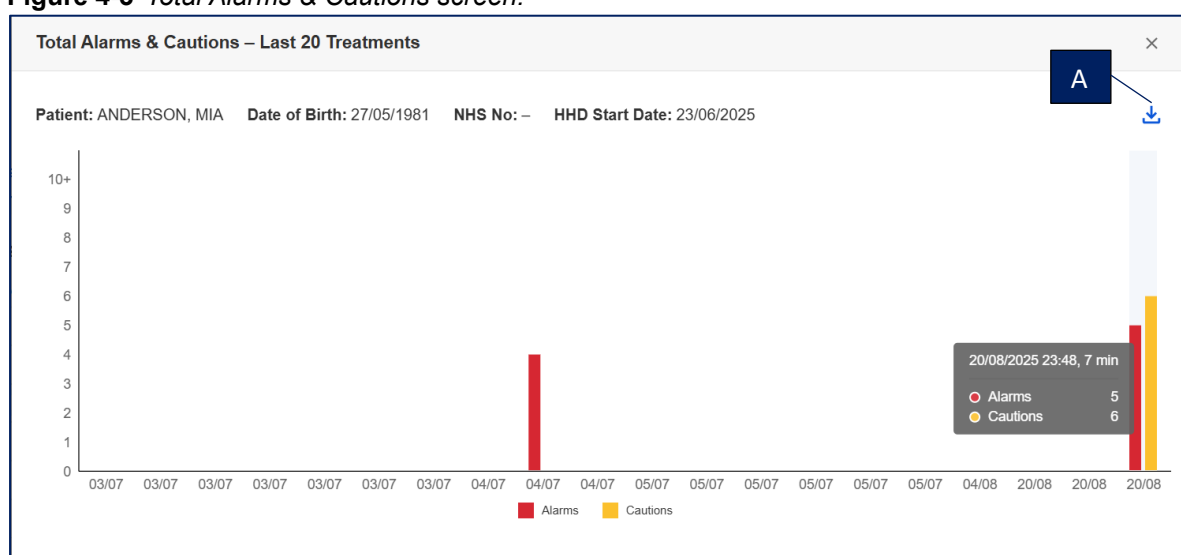

Hovering your mouse over the graph will display treatment date, duration, and total number of alarms and/or cautions.

Clicking on graph will take you to the Treatment Flowsheet screen. For more information, see [Chapter 5: Treatment Flowsheet](#).

Clicking on download icon, [Figure 4-3 \(A\)](#), allows you to download and save the Total Alarms & Cautions as a PDF, which can then be printed on a local printer.

Prescription Summary

The Prescription Summary provides an overview of the patient's therapy details, including prescribed therapy and supplies.

Clicking the → icon, [Figure 4-1 \(C\)](#), takes you to the *Prescription* section under the **Therapy Management** tab, where you can view and edit the patient's therapy details and supplies treatment plan. For more information, see the [Therapy Management](#) section.

Alert Criteria

The alert criteria entered and displayed in the patient profile determines the ranges at which a clinical alert (🚨) will be displayed.

While in the patient profile, the alert criteria card displays a preview of alert criteria thresholds that have been configured for blood pressure and temperature.

Clicking the → icon, [Figure 4-1 \(D\)](#), takes you to the *Alert Criteria* section under the **Therapy Management** tab, where you can view and edit the patient's blood pressure, temperature, and other alerts. For more information, see the [Therapy Management](#) section.

Lifestyle Goals

The Lifestyle goals card displays lifestyle goals for patients and/or care-partners as entered by clinicians into kinexus HHD. The date of the most recent update is displayed.

Clicking the ↗ icon, [Figure 4-1 \(E\)](#), opens a pop-up window, as shown in [Figure 4-4](#).

Once the patient and care partner details are entered, click the **Save** button, [Figure 4-4 \(A\)](#), to successfully save the information. Clicking **Cancel**, [Figure 4-4 \(B\)](#) or the **X**, [Figure 4-4 \(C\)](#), will take you back to the Patient Profile screen without making any changes to the lifestyle goals.

Figure 4-4 *Lifestyle Goals screen.*

Clicking on **See Latest Treatment**, [Figure 4-1 \(F\)](#), will take you to the Treatment Flowsheet. For more information, see [Chapter 5: Treatment Flowsheet](#).

Treatments

The **Treatments** tab, [Figure 4-5](#), lists treatments according to start date and time. Only treatments that your user role and designated clinic authorize you to view will be visible. To sort treatments from oldest to newest dates, click the downward arrow next to the **Date, Start** column header. For more details on the data displayed, see section, [Last 7 Days Treatments](#).


Figure 4-5 Patient profile – Treatments tab.

	Date, Start ↓	DUR, min	Alarms	Notes	Weight, kg	BP, mmHg	Pulse, bpm	UF, L	Dialysate, L	Blood, L
☐	20/08/2025, 23:48	7	6		62.3 pre 60 post 60 dry	120/80 pre 120/80 post	70 pre 70 post	1 actual 2.3 target -- rate	4 processed 2 prescribed	27
☐	20/08/2025, 23:41	0	1		62.4 pre -- post 60 dry	-- pre -- post	-- pre -- post	-- actual 2.4 target -- rate	0 processed 2 prescribed	
☐	20/08/2025, 23:38	0	6		100 pre 100 post 60 dry	120/80 pre 120/80 post	70 pre 70 post	-- actual 30 target -- rate	0 processed 2 prescribed	
☐	04/08/2025, 12:05	18			-- pre -- post 60 dry	-- pre -- post	-- pre -- post	0.1 actual -- target 0.46 rate	0.5 processed 30 prescribed	5.7
☐	05/07/2025, 12:31	0			-- pre -- post 50 dry	-- pre -- post	-- pre -- post	-- actual -- target -- rate	0 processed 40 prescribed	
☐	05/07/2025, 12:31	0			-- pre -- post 50 dry	-- pre -- post	-- pre -- post	-- actual -- target -- rate	0 processed 40 prescribed	
☐	05/07/2025, 12:31	0			-- pre -- post 50 dry	-- pre -- post	-- pre -- post	-- actual -- target -- rate	0 processed 40 prescribed	

Sort Treatment list

At the top of the *Treatments list* is a sort drop-down menu, [Figure 4-5](#) (A), which allows you to sort the list in four ways, as shown in [Figure 4-6](#).

Figure 4-6 Sort treatment drop-down list.

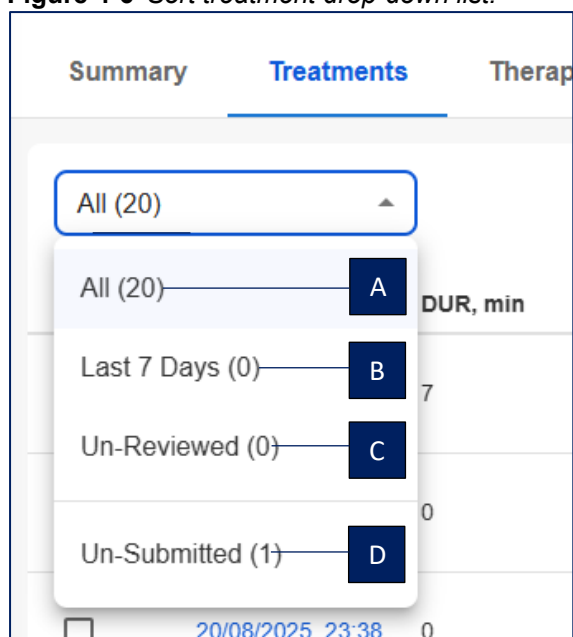

- **All**, [Figure 4-6 \(A\)](#), displays a complete list of treatments received for that patient.
- **Last 7 Days**, [Figure 4-6 \(B\)](#), displays treatments received within the last seven days.
- **Un-Reviewed**, [Figure 4-6 \(C\)](#), displays treatments that have not yet been reviewed.
- **Un-Submitted**, [Figure 4-6 \(D\)](#), displays treatments of that patient that have not yet been submitted and remain editable.

The treatment table, [Figure 4-5 \(B\)](#), displays a list of information related to each treatment. For more information on the treatment table, see the [Last 7 Days Treatments Table](#) section.

Clicking the checkbox next to a specific treatment date, enables the **Mark as Reviewed** or **Mark as Un-Reviewed**, [Figure 4-5 \(C\)](#), buttons, at the top of the table. For more information. See [How to Review Treatments](#) section.

Clicking on treatment will take you to the Treatment Flowsheet screen. For more information, see [Chapter 5: Treatment Flowsheet](#).

Therapy Management

The **Therapy Management** tab provides an overview of the patient's therapy details, including prescribed therapy, supplies, and configured alert criteria thresholds for blood pressure and temperature. This tab is divided into two sections:

- **Prescription**
- **Alert Criteria**

Prescription

The Prescription section provides an overview of the patient's therapy details, including prescribed therapy and supplies. To update information, click the field you want to edit.

After making the necessary changes in **Therapy** section, click the **Save** button, [Figure 4-7 \(A\)](#), to save the updated therapy details successfully. To discard edits and restore the original values, click **Clear**, [Figure 4-7 \(B\)](#).

Similarly, after making updates in the **Supplies** section, click the **Save** button, [Figure 4-7 \(C\)](#), to save the changes. To revert to the original values, click **Clear**, [Figure 4-7 \(D\)](#).

Figure 4-7 Therapy Management – Prescription screen.

The screenshot displays the 'Therapy Management' tab with the 'Prescription' sub-tab selected. The interface includes a top navigation bar with 'Summary', 'Treatments (• 11)', and 'Therapy Management'. The 'Prescription' section is active, showing 'Acknowledged by patient on 12/09/2025' and patient identifiers 'NxStage ID: 1073' and 'App Pin: 100319'. The 'Therapy' section contains input fields for 'Dry Weight *' (39.9 kg), 'Max. UFR *' (1.3 L/hr), 'Avg. UFR *' (1 L/hr), 'Initial Bloodflow Rate *' (560 ml/min), 'Target Bloodflow Rate *' (570 ml/min), and 'Dialysate Flow Rate *' (18). A 'Needle Gauge *' dropdown is set to 'None'. A 'Dosing Calculator' link is available. The 'Supplies' section includes 'Treatment Plan' (SAK-306(45Lact, 3Ca x 50L)), 'Medication' (AntiCoagulant), 'Modality' (Catheter), 'Dialysate Volume *' (56 L), 'Strength *' (17000), 'Treatments per Week *' (5 Days), 'Cartridge Type *' (Generic), and 'Route *' (Sub-Cutaneous). Both sections feature 'Clear' and 'Save' buttons. Update timestamps are shown at the bottom of each section: 'Updated by Birajan Shakya on 11/09/2025' for Therapy and 'Updated by Birajan Shakya on 12/09/2025' for Supplies. A note at the bottom states: 'After saving the therapy, the user's email ID and the date of the update will be displayed at the bottom of the screen.'



NOTE:

After saving the therapy, the user's email ID and the date of the update will be displayed at the bottom of the screen.

The **Therapy** section fields display the following information:

- **Dry Weight:** The target weight defined in the patient's prescription therapy in kilogrammes (kg).
- **Max. UFR:** The maximum ultrafiltration rate in litres per hour (L/hr).
- **Avg. UFR:** The average ultrafiltration rate in litres per hour (L/hr).
- **Initial Bloodflow Rate:** Initial (starting) blood flow rate prescribed for the patient, in millilitres per minute (ml/min).
- **Target Bloodflow rate:** Target blood flow rate prescribed for the patient, in millilitres per minute (ml/min).
- **Dialysate Flow Rate:** The rate of dialysate fluid flow in litres per hour (L/hr).
- **Needle Gauge:** The outer diameter (thickness) of the needle used for dialysis.

The **Supplies** section contains the following information/options:

- **Treatment Plan**
 - Solution Type
 - Dialysate Volume (in L)
 - Treatments per Week
 - Cartridge type
-
- **Medication**
 - Medication Type
 - Strength
 - Route
-
- **Modality**
 - Access Type
 - NxStage Type
-

The Supplies section of the prescription displays approved prescription values from Order Manager. Please refer to Order Manager documentation for more information on the prescription fields and process for prescription approval.

When entering values outside of the range set for that field, the field will be highlighted in red and save button will be disabled as shown in [Figure 4-8](#).

Figure 4-8 Incorrect value error screen.

Therapy [Dosing Calculator](#)

Dry Weight * 577 kg
Valid range: 0.1 - 300.0 (in increments of 0.1)

Max. UFR * 2.33 L/hr

Avg. UFR * 2.32 L/hr

Initial Bloodflow Rate * 600 ml/min

Target Bloodflow Rate * 550 ml/min

Dialysate Flow Rate * 16.7 L/hr

Needle Gauge * 16

Updated by Pavan Ramesha on 01/09/2025

[Clear](#) [Save](#)

Therapy Status

The Prescription Summary screen displays the status of the therapy in the top left corner. This feature provides clinicians with essential information regarding patient engagement, aiding in treatment monitoring and management.

- Initially, when clinicians add or modify the therapy, the status appears as **Not acknowledged by patient**, as shown in [Figure 4-9](#) (A).
- Once the patient has viewed and acknowledged the therapy change in the Nx2me App, the status changes to **Acknowledged by patient** along with the date, as shown in [Figure 4-10](#) (A).

Figure 4-9 Prescription Status – Not acknowledged by patient.

Patients (• 1) [Therapy Management](#)

Not acknowledged by patient **A**

Therapy

Dry Weight * 57 kg

Max. UFR * 2.33

Initial Bloodflow Rate * 600 ml/min

Target Bloodflow Rate * 550

Figure 4-10 Prescription Status - Acknowledged by patient.

Patients (• 1) [Therapy Management](#)

Acknowledged by patient on 01/09/2025 **A**

Therapy

Dry Weight * 49.1 kg

Max. UFR * 2.2

Initial Bloodflow Rate * 200 ml/min

Target Bloodflow Rate * 600

Alert Criteria

The Alert Criteria entered and displayed in the patient profile determines the ranges at which a clinical alert (🔔) will be displayed.

The Alert Criteria section provides an overview of the patient's alert criteria thresholds that have been configured for blood pressure and temperature.

Here clinicians can enter or edit the blood pressure, temperature, and other alerts for that individual patient. After making the necessary changes, click the **Save** button, [Figure 4-11 \(A\)](#), to save the updated therapy details successfully. To discard edits and restore the original values, click **Clear**, [Figure 4-11 \(B\)](#).

Figure 4-11 Therapy Management – Alert Criteria screen.

When entering values outside of the range set for that field, the field will be highlighted in red and save button will be disabled, as shown in [Figure 4-12](#).

Figure 4-12 Incorrect value error screen.



NOTE:

After saving the alert criteria, the user's email ID and the date of the update will be displayed at the bottom of the screen.

The following alert criteria can be defined:

Blood Pressure

- **Min. Systolic:** The minimum systolic blood pressure threshold, in mmHg, allowed before an alert is displayed.
- **Max. Systolic:** The maximum systolic blood pressure threshold, in mmHg, allowed before an alert is displayed.
- **Min. Diastolic:** The minimum diastolic blood pressure threshold, in mmHg, allowed before an alert is displayed.
- **Max. Diastolic:** The maximum diastolic blood pressure threshold, in mmHg, allowed before an alert is displayed.

Temperature

- **Max:** Maximum temperature threshold in the temperature units specified for the patient.

Other Alerts

- **Dialysate Vol. Deviation:** The volume of dialysate used in a treatment deviate from the prescribed volume by more than the specified amount in litres (L).
- **Dry Weight Deviation:** The patient's dry weight at the end of a treatment deviates from the prescribed dry weight by more than the specified amount in kilogrammes (kg).
- **No. of Cycler Alarms:** The total number of cycler alarms is greater than or equal to the specified value.



NOTE:

The alert criteria defined for Blood Pressure and Temperature will be displayed in the Nx2me App as an alert if a value is that exceeds the prescribed range. Please see the Nx2me App user guide for more information.


END OF CHAPTER

Chapter 5: Treatment Flowsheet

Overview

A treatment flowsheet is an electronic record sent from the patient app (Nx2me) to the Provider Portal that provides details about each dialysis treatment. Each flowsheet is a combination of data collected directly from the cyclor and data entered by the patient into the Nx2me app.

Figure 5-1 *Flowsheet Overview.*

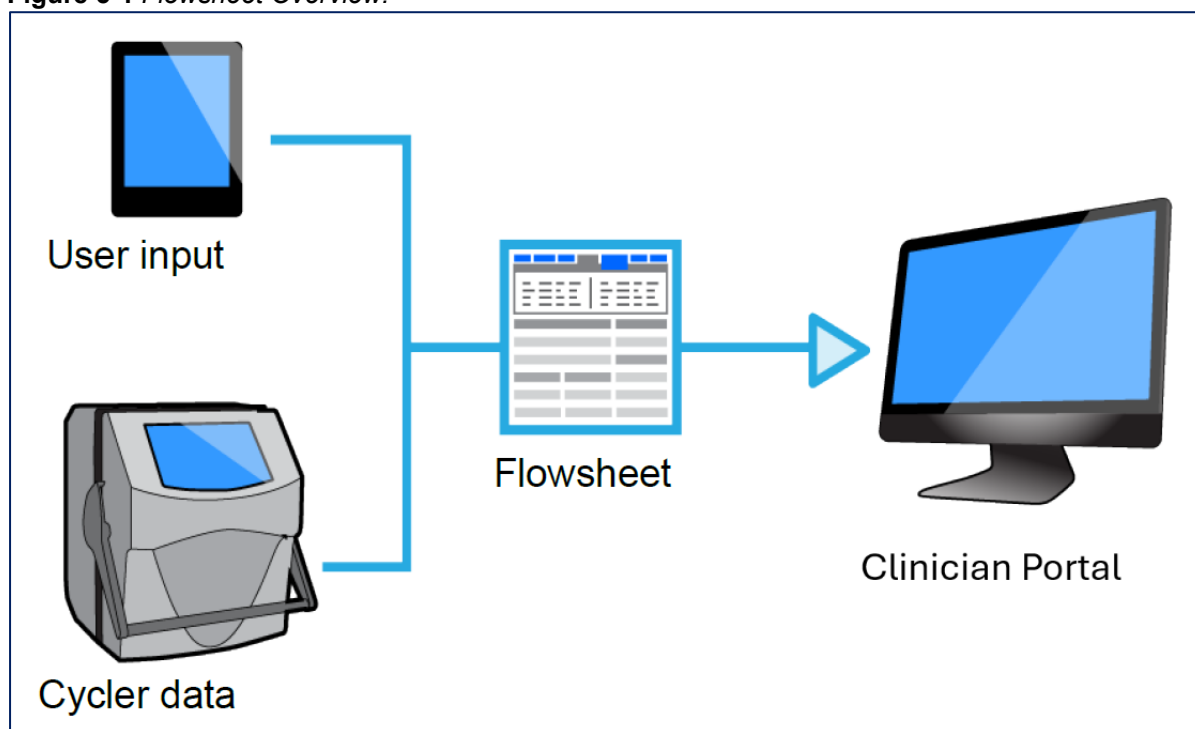

Treatment Status

The Flowsheet screen displays the status of the treatment flowsheet in the top left corner of the screen.

- Once the patient confirms the treatment in the Nx2me app, the status appears as **Submitted by patient** along with the date and timestamp, as shown in [Figure 5-2](#) (A).
- If the patient has not confirmed the treatment that status will be displayed as **Un-Submitted**, as shown in [Figure 5-3](#) (A).
 - When flowsheets are in an un-submitted state they are still editable by the patient
 - Flowsheets that have a dialysate volume processed and treatment duration associated with them will automatically be submitted after 3 days. Flowsheets that have been automatically submitted will display **Submitted by System** along with the date and timestamp.



NOTE:

Internet connectivity is required for flowsheet submission to the cloud and is required for both confirmation by the patient and system confirmation.

Figure 5-2 Flowsheet Status – Submitted by patient.

Clinic End to End 2 / AJOHN, ADAM / Treatment Flowsheet

01/09/2025 – 0 0 0

Start Time: 15:05 | Duration: 0 min | Submitted by patient: 01/09/2025, 15:06 **A**

Summary | Cyclor Data | Vitals, Medication & Equipment

Figure 5-3 Flowsheet Status – Un-Submitted.

Clinic End to End 2 / DemoPatient D, Steve / Treatment Flowsheet

03/09/2025 – 0 0 0

Start Time: 20:05 | Duration: 0 min | Un-Submitted **A**

How to access the treatment flowsheet

Clicking on a treatment from the list will open its corresponding Treatment Flowsheet. Here, you can view the flowsheet details, such as start time, duration, and status.

The Treatment Flowsheet screen includes the following tabs to easily access its main functions.

- **Summary** – Provides an overview of the treatment; see [Summary](#) for more information.
- **Cycler Data** – Displays cycler-related data; see [Cycler Data](#) for more information.
- **Vitals, Medication & Equipment** – Displays vital signs, medications, and equipment information. See [Vitals, Medication & Equipment](#) for more information.

Summary

The **Treatment Flowsheet – Summary** tab contains the following sections:

- **Dialysate**, [Figure 5-4](#) (A)
- **Ultrafiltration**, [Figure 5-4](#) (B)
- **Blood**, [Figure 5-4](#) (C)
- **Alerts, Alarms & Cautions**, [Figure 5-4](#) (D)
- **Weight**, [Figure 5-4](#) (E)
- **Prescription**, [Figure 5-4](#) (F)
- **Clinic Note**, [Figure 5-4](#) (G)
- **Vitals**, [Figure 5-4](#) (H)
- **Medication**, [Figure 5-4](#) (I)
- **Patient Notes & Photos**, [Figure 5-4](#) (J)

Figure 5-4 *Treatment Flowsheet screen.*

PRESENIUS MEDICAL CARE

Clinic End to End 2 / ANDERSON, MIA / Treatment Flowsheet

01/09/2025

4 9

Start Time: 08:59

Duration: 40 min

Submitted by patient: 01/09/2025, 09:43

Mark as Un-Reviewed

Download

Summary

Cycler Data

Vitals, Medication & Equipment

Dialysate

23%

Processed, L

1.2

Prescribed, L

5.2

Ultrafiltration

Target, L

3.1

Actual, L

1.2

Average, L/hr

1.55

Blood

Processed, L

10.6

Alerts, Alarms & Cautions

9 Alerts

4 Warnings

4 Errors

Clinical Alerts

Pre Temperature, °C

35.7

Post Temperature, °C

35.9

Post-Dialysis Weight, kg

60.2

Number of cycler alarms exceeded Alert Criteria threshold

-

Dialysate Processed, L

1.2

During treatment - BP, mmHg

01/09/2025 09:11 - 121/120

During treatment - BP, mmHg

01/09/2025 09:12 - 120/121

During treatment - BP, mmHg

01/09/2025 09:13 - 69/70

During treatment - BP, mmHg

01/09/2025 09:13 - 70/69

Cycler Alarms

A11

1

A10

1

Weight

Pre-Dialysis Weight, kg

61.1

Weight to Remove, kg

3.1

Post-Dialysis Weight, kg

60.2

Dry Weight (Rx), kg

58

Prev Post-Weight, kg

60.1

Prescription

Prescription was changed on 01/09/2025

Dry Weight, kg

58

Max. UFR, L/hr

2.33

Avg. UFR, L/hr

2.32

Dialysate Flow Rate, L/hr

16.7

Initial Blood Flow Rate, ml/min

600

Target Blood Flow Rate, ml/min

550

Needle Gauge

16

Clinic Note

Enter your notes here for reference.

Updated by Sandhya Kotturi on 01/09/2025

Vitals 2

Medication 5

Patient Notes & Photos

Pre-Treatment

09:43

The European languages are members of the same family. Their separate existence is a myth. For science, music, sport, etc, Europe uses the same vocabulary. The languages only differ in their grammar.

3 Photos

Treatment

09:13

The European languages are members of the same family. Their separate existence is a myth. For science, music, sport, etc, Europe uses the same vocabulary. The languages only differ in their grammar.

0 Photos

Post-Treatment

09:43

But I must explain to you how all this mistaken idea of denouncing pleasure and praising pain was born and I will give you a complete account of the system, and expound the actual teachings of the great explorer of the truth, the master-builder of human happiness. No one rejects, dislikes, or avoids pleasure itself, because it is pleasure, but because those who do not know how to pursue pleasure rationally encounter consequences that are extremely painful. Nor again is there anyone who loves or


3 Photos

Dialysate

The Dialysate section displays the amount of dialysate fluid processed during the treatment and the amount prescribed by the patient's physician in litres (L), as shown in [Figure 5-5](#). It includes the following measurements:

- **Processed, L:** The cyclor dialysate volume processed during the treatment, in litres (L).
- **Prescribed, L:** The volume of dialysis prescribed to use for each treatment, in litres (L).

Figure 5-5 *Dialysate.*

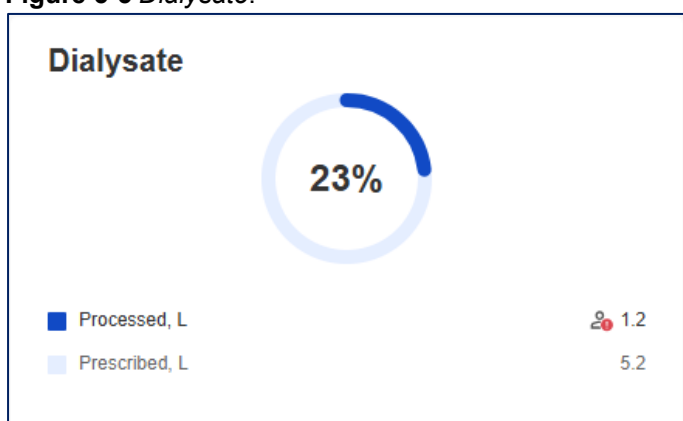

Ultrafiltration

The Ultrafiltration section displays the target and actual litres (L) of ultrafiltration removed from the patient during treatment as well as the average ultrafiltration rate in litres per hour (L/hr), as shown in [Figure 5-6](#). It includes the following measurements:

- **Target, L:** The prescribed ultrafiltration volume to be removed during the treatment, in litres (L).
- **Actual, L:** The ultrafiltration volume removal achieved during the treatment in litres (L).
- **Average, L/hr:** The average rate of ultrafiltration, in litres per hour (L/hr).


Figure 5-6 *Ultrafiltration.*

Ultrafiltration	
Target, L	3.1
Actual, L	1.2
Average, L/hr	1.55

Blood

The Blood section displays the total amount of blood processed in litres (L), as shown in [Figure 5-7](#).

Figure 5-7 *Blood.*

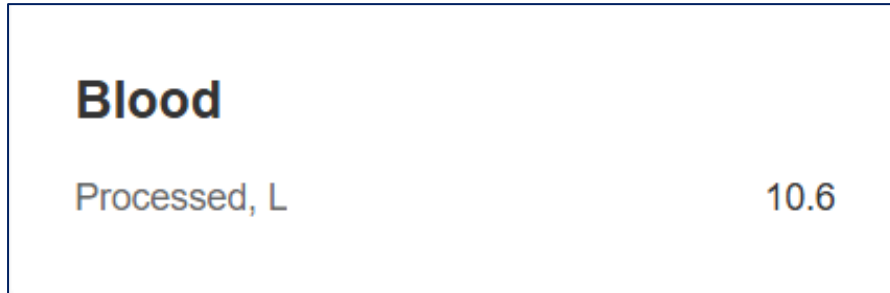

Alert, Alarms & Cautions

The **Alerts, Alarms, and Cautions** card displays clinical alerts, cyclor alarms, and cyclor cautions recorded for an individual treatment, as shown in [Figure 5-8](#). It includes the following measurements:

- **Pre Temperature, °C** – Patient's body temperature before starting the treatment, in Celsius (°C).
- **Post Temperature, °C** - Patient's body temperature after completing the treatment, in Celsius (°C).
- **Post-Dialysis Weight, kg** - Patient's weight after completing the treatment, in kilograms (kg).
- **Number of cyclor alarms exceeded Alert Criteria threshold** – The total number of cyclor alarms that exceeded the alert criteria threshold.
- **Dialysate Processed, L** – Displays the cyclor dialysate volume processed during the treatment, in litres (L).
- **During treatment – BP, mmHg** - Patient's blood pressure recorded during the treatment, in millimeters of mercury (mmHg).
- **Cyclor Alarms** - List of cyclor alarm codes triggered during the treatment.
- **Cyclor Cautions** - List of cyclor caution codes triggered during the treatment.

Figure 5-8 Alert, Alarms & Cautions.

Alerts, Alarms & Cautions		🔔 9	⚠️ 4	🟡 4
🔔 Clinical Alerts				
Pre Temperature, °C				35.7
Post Temperature, °C				35.9
Post-Dialysis Weight, kg				60.2
Number of cyclor alarms exceeded Alert Criteria threshold				—
Dialysate Processed, L				1.2
During treatment - BP, mmHg	01/09/2025 09:11 - 121/120			
During treatment - BP, mmHg	01/09/2025 09:12 - 120/121			
During treatment - BP, mmHg	01/09/2025 09:13 - 69/70			
During treatment - BP, mmHg	01/09/2025 09:13 - 70/69			
⚠️ Cyclor Alarms				🔔
A11				1
A10				1

Weight

The Weight section displays weight data entered by the patient before and after treatment, as well as the weights used to calculate treatment targets, as shown in [Figure 5-9](#).

- **Pre-Dialysis Weight, kg:** The weight of the patient before starting the treatment, in kilogrammes (kg).
- **Weight to Remove, kg:** The calculated amount of weight to be removed from the patient during treatment in kilogrammes (kg).
- **Post-Dialysis Weight, kg:** The weight of the patient after completing the treatment, in kilogrammes (kg).
- **Dry Weight (Rx), kg:** The patient's prescribed target weight post-dialysis, in kilogrammes (kg).
- **Prev. Post-Weight, kg:** The patient's previous post-dialysis weight. If available from previous flowsheets, this value is captured by the Nx2me app and added to the flowsheet automatically; otherwise, it is entered by the patient, in kilogrammes (kg).

Figure 5-9 Weight.

Weight	
Pre-Dialysis Weight, kg	61.1
Weight to Remove, kg	3.1
Post-Dialysis Weight, kg	🔔 60.2
Dry Weight (Rx), kg	58
Prev. Post-Weight, kg	60.1

Prescription

The Prescription section displays information about the patient's prescription settings on the date of the treatment, as shown in [Figure 5-10](#). It includes the following details:

- **Dry Weight, kg:** The target weight in kilogrammes (kg).
- **Max. UFR, L/hr:** The maximum ultrafiltration rate in litres per hour (L/hr).
- **Avg. UFR, L/hr:** The average ultrafiltration rate in litres per hour (L/hr).
- **Dialysate Flow Rate, L/hr:** The rate of dialysate fluid flow in litres per hour (L/hr).
- **Initial Blood Flow Rate, ml/min:** The initial (starting) rate of blood flow in millilitres per minute (ml/min).
- **Target Blood Flow Rate, ml/min:** The target rate of blood flow in millilitres per minute (ml/min).
- **Needle Gauge:** The outer diameter (thickness) of the needle used for dialysis.

Figure 5-10 Therapy Prescription.

Prescription	
Prescription was changed on 01/09/2025	
Dry Weight, kg	58
Max. UFR, L/hr	2.33
Avg. UFR, L/hr	2.32
Dialysate Flow Rate, L/hr	16.7
Initial Blood Flow Rate, ml/min	600
Target Blood Flow Rate, ml/min	550
Needle Gauge	16

Clinic Note

The Clinic Note card displays clinical notes entered by clinicians into kinexus HHD. Clicking the  icon, [Figure 4-1](#) (E), opens a pop-up window, as shown in [Figure 5-11](#).

Once the clinic note details are entered, click the **Save** button, [Figure 5-11](#) (A), to successfully save the information. Clicking **Cancel**, [Figure 5-11](#) (B) to discard edits and restore the original notes. Clicking the **X**, [Figure 5-11](#) (C), will take you back to the *Treatment Flowsheet* screen without saving any changes.

Figure 5-11 Clinic Note.

**NOTE:**

After saving the clinic notes, the user's email ID and the date of the update will be displayed at the bottom of the screen.

Vitals

The **Vitals** section displays vital sign data entered by the patient. Clicking the  icon, [Figure 5-4 \(H\)](#), takes you to the **Vitals, Medication & Equipment** tab, where you can view detailed vital sign information. For more information, see the [Vitals, Medication & Equipment](#) section.

Medication

The **Medication** section displays medications administered during treatment. Clicking the  icon, [Figure 5-4 \(I\)](#), takes you to the **Vitals, Medication & Equipment** tab, where you can view detailed medication information. For more information, see the [Vitals, Medication & Equipment](#) section.

Patient Notes & Photos

The Patient Notes & Photos section displays notes or images entered by the patient for this treatment, as shown in [Figure 5-12](#). This section includes Pre-Treatment, Treatment, and Post-Treatment notes and photos.

- **Pre-Treatment:** Notes and photos entered before the treatment.
- **Treatment:** Notes and photos were entered during the treatment.
- **Post-Treatment:** Notes and photos entered after the treatment.

Figure 5-12 Patient Notes & Photos.

Patient Notes & Photos	Treatment	Post-Treatment
Pre-Treatment  09:43 The European languages are members of the same family. Their separate existence is a myth. For science, music, sport, etc, Europe uses the same vocabulary. The languages only differ in their grammar.  3   	Treatment  09:13 The European languages are members of the same family. Their separate existence is a myth. For science, music, sport, etc, Europe uses the same vocabulary. The languages only differ in their grammar.  0 No photos  09:14 But I must explain to you how all this mistaken idea of denouncing pleasure and praising pain was born and I will give you a complete account of the system, and expound the actual teachings of the great explorer of the truth, the master-builder of human happiness. No one rejects, dislikes, or avoids pleasure itself, because it is pleasure, but because those who do not know how to pursue pleasure rationally encounter consequences that are extremely painful. Nor again is there anyone who loves or.  3   	Post-Treatment  09:43 But I must explain to you how all this mistaken idea of denouncing pleasure and praising pain was born and I will give you a complete account of the system, and expound the actual teachings of the great explorer of the truth, the master-builder of human happiness. No one rejects, dislikes, or avoids pleasure itself, because it is pleasure, but because those who do not know how to pursue pleasure rationally encounter consequences that are extremely painful. Nor again is there anyone who loves or.  3   

Clicking on the image opens a pop-up window with an enlarged view, as shown in [Figure 5-13](#).

Figure 5-13 Patient Photos.



Review and Un-Review Treatment

The **Mark as Un-Reviewed** button, [Figure 5-14](#) (A), appears if you select a treatment that has been reviewed. When a treatment flowsheet is opened, the status automatically changes to reviewed. Users can review or unreview treatment on this screen. For more information, see the [How to Review Treatments](#) section in [Chapter 3: Clinic Dashboard](#).

Figure 5-14 Mark as Un-Reviewed.

Clinic End to End 2 / ANDERSON, MIA / Treatment Flowsheet		B	
01/09/2025 –  4  9  1		← 20/08/2025	
Start Time: 08:59 Duration: 40 min Submitted by patient: 01/09/2025, 09:43		A  	
Summary Cycler Data Vitals, Medication & Equipment			

Clicking on download icon, [Figure 5-14](#) (B), allows you to download and save the Treatment Flowsheet as a PDF, which can then be printed on a local printer.

Cycler Data

The **Treatment Flowsheet – Cycler Data** tab displays all the cycler data recorded during Prime, Treatment, and Rinseback, as shown in [Figure 5-15](#). It also includes alarms, cautions, and vitals if any were recorded or entered during the treatment.

- **Alarms:** Displays the alarms that occurred. Alarm data is only included in the flowsheet if the patient is running the Nx2me app, [Figure 5-15 \(A\)](#).
- **Cautions:** Displays the cautions that occurred, as shown in [Figure 5-15 \(B\)](#).
- **Vitals:** Displays the vital signs entered by the patient, [Figure 5-15 \(C\)](#).

Figure 5-15 Treatment table.

Summary Cycler Data Vitals, Medication & Equipment								
Number of cycler alarms exceeded Alert Criteria threshold								
Dialysate			Ultrafiltration					
Time	Rate, L/hr	Vol, L	Rate, L/hr	Vol, L	BFR, mL/min	AP, mmHg	VP, mmHg	EP, mmHg
Prime								
06:41	A11							
06:48	C80							
06:48	C81							
Treatment								
09:01	1.8	0.1	2.32	0	260	—	100	230
09:11	Blood pressure, mmHg: 120/120 Pulse, bpm: 70 Access Status: Good Saline Bolus, mL: 50							
09:11	Blood pressure, mmHg: 121/120 Pulse, bpm: 70 Access Status: Good Saline Bolus, mL: 50							
09:12	Blood pressure, mmHg: 120/121 Pulse, bpm: 70 Access Status: Good Saline Bolus, mL: 50							
09:12	Blood pressure, mmHg: 70/70 Pulse, bpm: 70 Access Status: Good Saline Bolus, mL: 50							
09:13	Blood pressure, mmHg: 69/70 Pulse, bpm: 70 Access Status: Inflamed Saline Bolus, mL: 50							
09:13	Blood pressure, mmHg: 70/69 Pulse, bpm: 70 Access Status: Swelling Saline Bolus, mL: 100							
09:23	C4							
09:29	1.8	0.9	2.32	0.9	260	—	100	230
09:38	A11							
09:39	A10							
09:41	0	1.2	0	1.2	260	—	100	230
Rinseback								
09:42	C4							
09:42	A11							

Cycler Data Table

The following information is displayed in the *Cycler Data table*:

- **Time:** The time the cycler data values or vitals were recorded. All data in the row was recorded at the indicated time.
- **Dialysate Rate, L/hr:** Cycler dialysate flow rate, in litres per hour (L/hr).
- **Dialysate Vol, L:** Cycler dialysate volume, in litres (L), that has been processed.
- **Ultrafiltration Rate L/hr:** Cycler ultrafiltration flow rate, in litres per hour (L/hr).
- **Ultrafiltration Vol, L:** Cycler ultrafiltration volume, in litres (L), that has been removed from the patient.
- **BFR, mL/min:** Cycler blood flow rate, in millilitres per minute (ml/min).
- **AP, mmHg:** Cycler arterial pressure, in millimetres of mercury (mmHg).
- **VP, mmHg:** Cycler venous pressure, in millimetres of mercury (mmHg).
- **EP, mmHg:** Cycler effluent pressure, in millimetres of mercury (mmHg).
- **Blood pressure:** The patient's blood pressure, in the format systolic/diastolic, in millimetres of mercury (mmHg), at the time indicated.
- **Pulse:** The patient's pulse at the time indicated, in beats per minute (bpm).
- **Access condition:** The patient's reported condition of their vascular access:
 - **Good**
 - **Inflamed**
 - **Swelling**
- **Saline Bolus:** The volume of the saline bolus, in millilitres (ml), that the patient administered at the time indicated. If the value of this field is 0, no saline bolus was administered at the time indicated.

Vitals, Medication & Equipment

The **Treatment Flowsheet – Vitals, Medication & Equipment** tab displays detailed information about the patient's vital signs, medications, and equipment data for that treatment.

Vitals

The Vitals section displays vital sign data entered by the patient, as shown in [Figure 5-16](#). Vital signs include Blood pressure (mmHg), Pulse (bpm), and Temperature (°C). This section only displays pre-treatment and post-treatment data.

- **Pre-Treatment:** Vital sign data entered by the patient before the treatment.
- **Post-Treatment:** Vital sign data entered by the patient after the end of the treatment.

Figure 5-16 Vitals.

Vitals		
Vital	Pre-Treatment	Post-Treatment
Blood pressure, mmHg	100/74	120/70
Pulse, bpm	66	70
Temperature, °C	35.7	35.9

Medication

The Medication section displays medications administered during the treatment, as entered by the patient in the Nx2me app, [Figure 5-17](#). It includes the following details:

- **Time Administered:** Time the medication was administered, as entered by the patient in the Nx2me app.
- **Name:** The name of the medication administered, as entered by the patient in the Nx2me app.
- **Dose:** The dose of the medication administered, as entered by the patient in the Nx2me app.
- **Route:** The location where the medication was administered, as entered by the patient in the Nx2me app.

Figure 5-17 Medication.

Medication			
Time Administered	Name	Dose	Route
09:08	Anticoagulant	134	Via Cartridge
09:08	ESA	36	Via Access
09:08	ESA	Anticoagulant	Oral
09:08	Clinic	074	Subcutaneous
09:15	Anticoagulant	69	Via Access

Equipment Data

The Equipment Data section displays information about the equipment and disposables used in the treatment, as entered by the patient, shown in [Figure 5-18](#). It includes Fluid Source, Cartridge Type, Cartridge Lot #, PAK Lot #, SAK Lot #, SAK Use #, and Chloramines Passed.

- **Fluid Source:** The therapy fluid source used during this treatment, either **PureFlow** or **bags**.
- **Cartridge Type:** The type of cartridge used during the treatment.
- **Cartridge Lot #:** The lot number of the cartridge used during the treatment.

The following fields are only displayed if the fluid source selected is **PureFlow**.

- **PAK Lot #:** Lot number of the PAK used during the treatment.
- **SAK Lot #:** Lot number of the SAK used during the treatment.
- **SAK Use #:** The number of times the SAK has been used for treatment, including this session.
- **Chloramines Passed?:** Confirmation that the chloramine level was tested prior to a treatment when a new SAK batch is made.

The following field is only displayed if the fluid source selected is **Bags**.

- **Warmer:** Indicates whether a fluid warmer was used during the treatment.

Figure 5-18 *Equipment Data.*

Equipment Data			
Fluid Source	PureFlow™ SL	SAK Lot #	07
Cartridge Type	CAR 125-B	SAK Use #	3
Cartridge Lot #	09	Chloramines Passed?	No
PAK Lot #	08		

END OF CHAPTER

Chapter 6: User Management

Overview

User management tasks can be performed by users with an admin role. Within the user management section, Admin users may create new users, assign user roles, or view and edit existing user roles and information.

Once the Menu icon, [Figure 6-1 \(A\)](#), is clicked from the side navigation bar, the screen-overlay menu appears, as shown in [Figure 6-2](#). Clicking **Users** button, [Figure 6-2 \(A\)](#), takes the user to the User Management page where a list of users is displayed. Clicking **kinexus HHD** button, [Figure 6-2 \(B\)](#), will take user back to the kinexus HHD portal.

Figure 6-1 Menu screen.

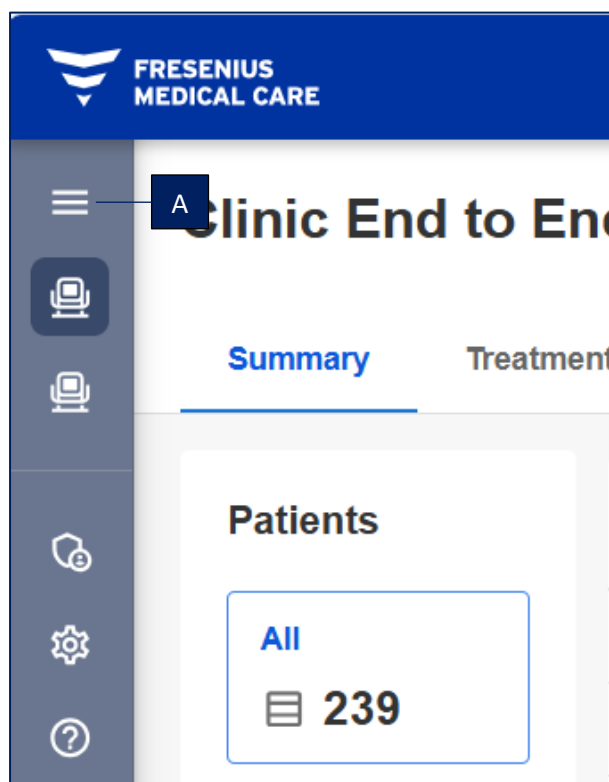
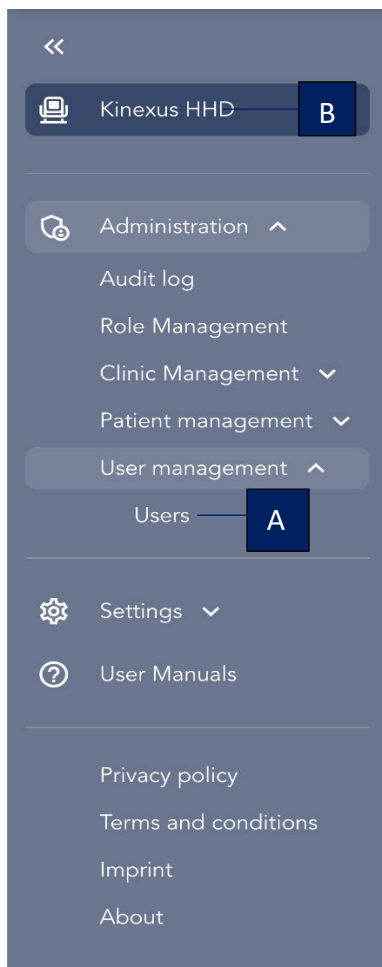


Figure 6-2 User management drop-down screen.

User Information

The Users table, [Figure 6-3](#) (A), provides the following user account information:

- **Name** displays the user's first name and last name.
- **Status** specifies whether the user is New, Draft, Active, or Inactive.
 - **New** indicates a user account was created, but the setup has not been completed by the user.
 - **Draft** indicates that particular user role has been assigned by the admin which activates the user, but the user has never logged in for the first time.
 - **Active** indicates the user has logged into the application, accepted the legal documents, and started using the application.
 - **Inactive** user accounts are inactive either manually or due to inactivity by the user (e.g., have not logged into the kinexus HHD portal in the last 3 months).
 - **Locked** user status implies that there has been a security measure applied to the account, due to multiple failed login attempts. Locked accounts may require special actions (contacting support or waiting for a specific period) to unlock and regain access. If user tries to log in during this state, an error message will appear that the account is locked.

- **Email** displays the user's email address.
- **Role** displays the role they are assigned within their organisation.

Figure 6-3 *User list screen.*



At the top of the user screen, there are four tabs: All, Active, Inactive, and Locked.

- **All:** This tab displays all users, regardless of their status.
- **Active:** This tab displays only users with an active status.
- **Inactive:** This tab displays only users with an inactive status.
- **Locked:** This tab displays only users with a locked status.

User Search

The search bar [Figure 6-3](#) (B). on top of user list allows you to search for users by first, last name, or email address. Searching happens automatically once you enter the letters in the search bar. To return to the full user list, remove your search parameters from the search bar.

Create User

Clicking **Create User** from the user list screen, [Figure 6-3](#) (C), will open a create user screen, [Figure 6-4](#). Here, you can enter the new user details, such as First name, Last name, Email Address, Phone Prefix (i.e., Country Code) with Phone Number, Language and select the preferred clinic from the drop-down list. Additionally, you can assign the necessary roles to the user under the Roles. **Kinexus Home HD User** role must be selected for each user to access the kinexus HHD application.

Figure 6-4 Create User screen.


Upon entering all the user details and clicking **Create**, [Figure 6-4](#) (A), it successfully saves the user information and sends an email to the user's email address with instructions to login to the kinexus HHD application. Clicking **Cancel**, [Figure 6-4](#) (B), opens a pop-up window as shown in [Figure 6-5](#).



NOTE:

The following fields in the Create User screen are mandatory for the user to be saved: First name, Last name, Email and Language. Additionally, assignment of roles to the user is necessary for the user to receive email with information about login and account set up.

Figure 6-5 Discard changes screen.

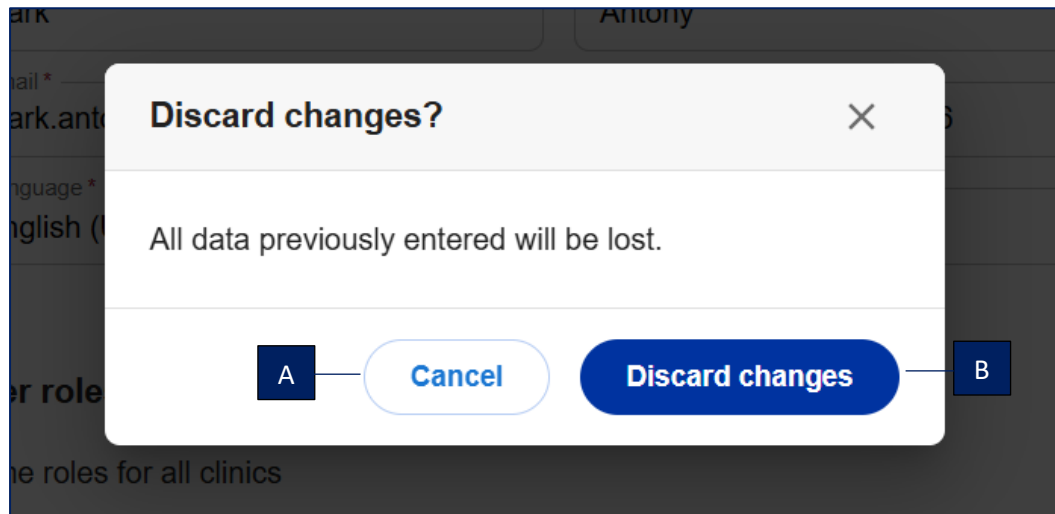

Clicking **Cancel**, [Figure 6-5](#) (A), takes you back to the *Create User* screen and clicking **Discard changes**, [Figure 6-5](#) (B), take you back to the users list screen.

Description of User Roles


kinexus HHD includes the following user account types (roles):

- **User Administrator**
 - Can create HHD user and assign the roles and clinics to the user.
 - Can read/modify data from HHD UI.
 - Can access all clinics as granted access by the system administrator (customer service support required).
- **Kinexus Home HD User**
 - Can read/modify data from HHD UI.
 - Can only access clinics as granted access by the Customer Administrator.
 - Can access both PD and HHD applications.
- **Read only User**
 - Can only view Patient Data and Prescription.
 - Can only view Treatments and Flowsheet details.
- **Patient Administartor**
 - Can add patients to kinexus HHD application.

Edit User Information

For any user, clicking on the ellipsis (⋮) symbol adjacent to their name in the user table list opens a screen-overlay menu, as shown in Figure 6-6. Clicking on **Edit**, Figure 6-6 (A), takes you to the edit user screen (Figure 6-7).

Figure 6-6 All Users screen.

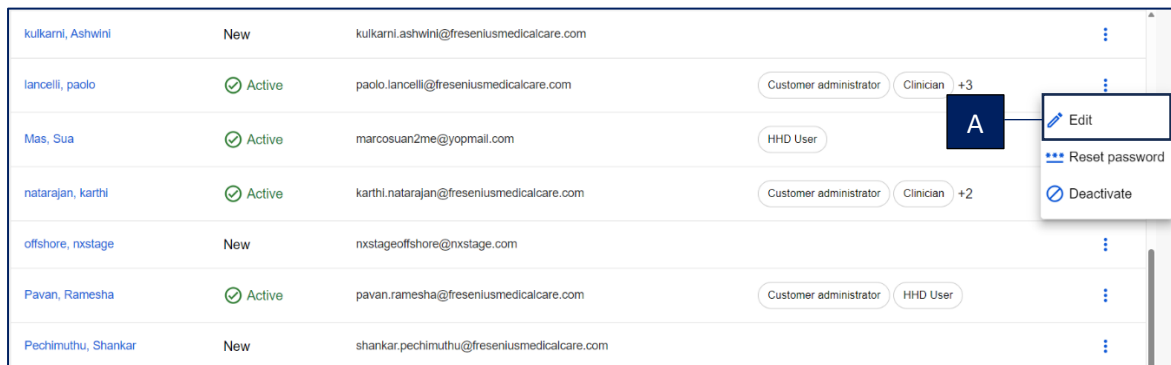


kulkarni, Ashwini	New	kulkarni.ashwini@freseniusmedicalcare.com		
lancelli, paolo	Active	paolo.lancelli@freseniusmedicalcare.com	Customer administrator	Clinician +3
Mas, Sua	Active	marcosuan2me@yopmail.com	HHD User	
natarajan, karthi	Active	karthi.natarajan@freseniusmedicalcare.com	Customer administrator	Clinician +2
offshore, nxstage	New	nxstageoffshore@nxstage.com		
Pavan, Ramesha	Active	pavan.ramesha@freseniusmedicalcare.com	Customer administrator	HHD User
Pechimuthu, Shankar	New	shankar.pechimuthu@freseniusmedicalcare.com		

Ellipsis menu for 'Mas, Sua':

- Edit (A)
- Reset password
- Deactivate

Figure 6-7 Edit User screen.

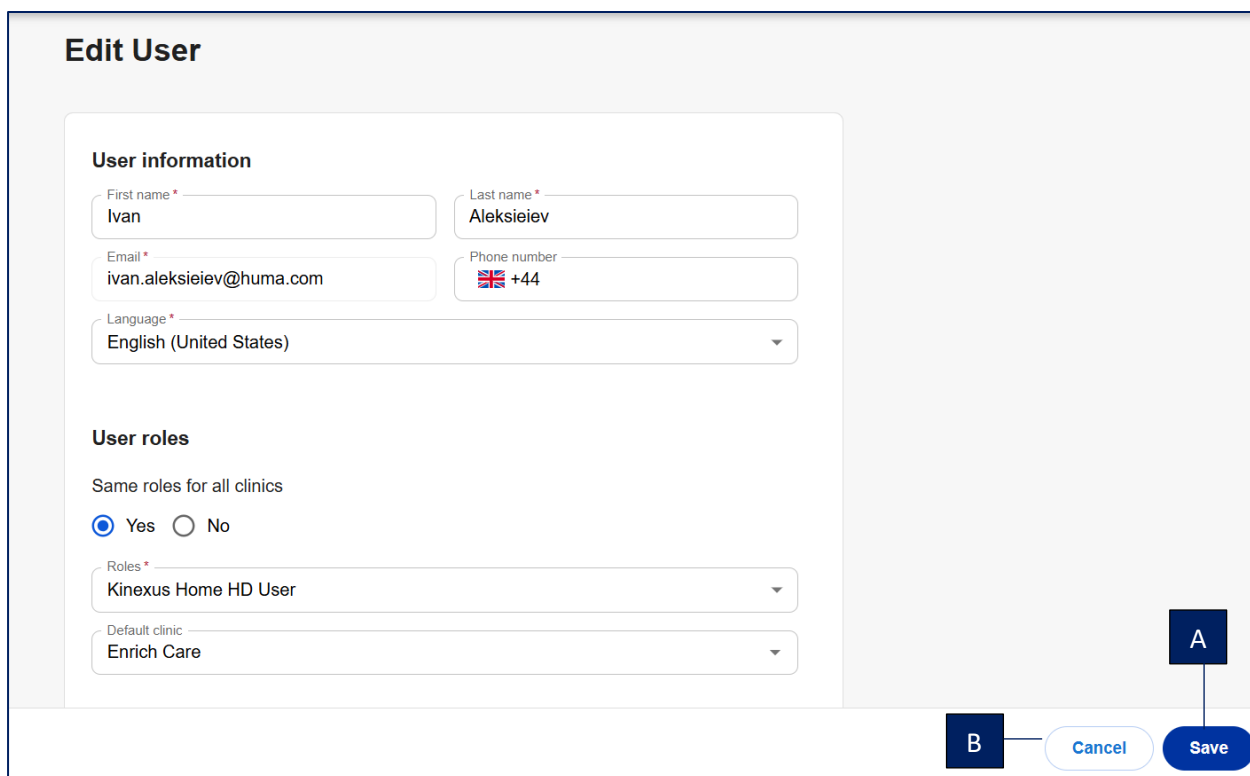

Edit User

User information

First name *
Ivan

Last name *
Aleksieiev

Email *
ivan.aleksieiev@huma.com

Phone number
+44

Language *
English (United States)

User roles

Same roles for all clinics

☒ Yes ☐ No

Roles *
Kinexus Home HD User


Default clinic
Enrich Care

A

B Cancel Save

Edit user screen allow you to edit the information, default clinic and roles of the user. Upon editing the user details and clicking **Save**, Figure 6-7 (A), it successfully saves the user information. Clicking **Cancel**, Figure 6-7 (B), opens a pop-up window as shown in Figure 6-8.

Figure 6-8 Discard changes screen.

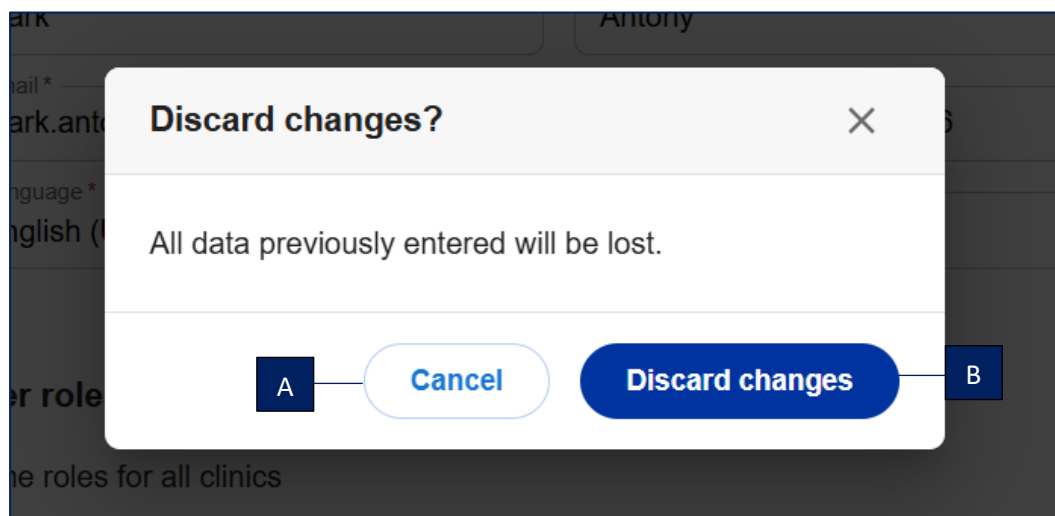

Clicking **Cancel**, [Figure 6-8 \(A\)](#), takes you back to the *Edit User* screen and clicking **Discard changes**, [Figure 6-8 \(B\)](#), take you back to the users list screen.

END OF CHAPTER